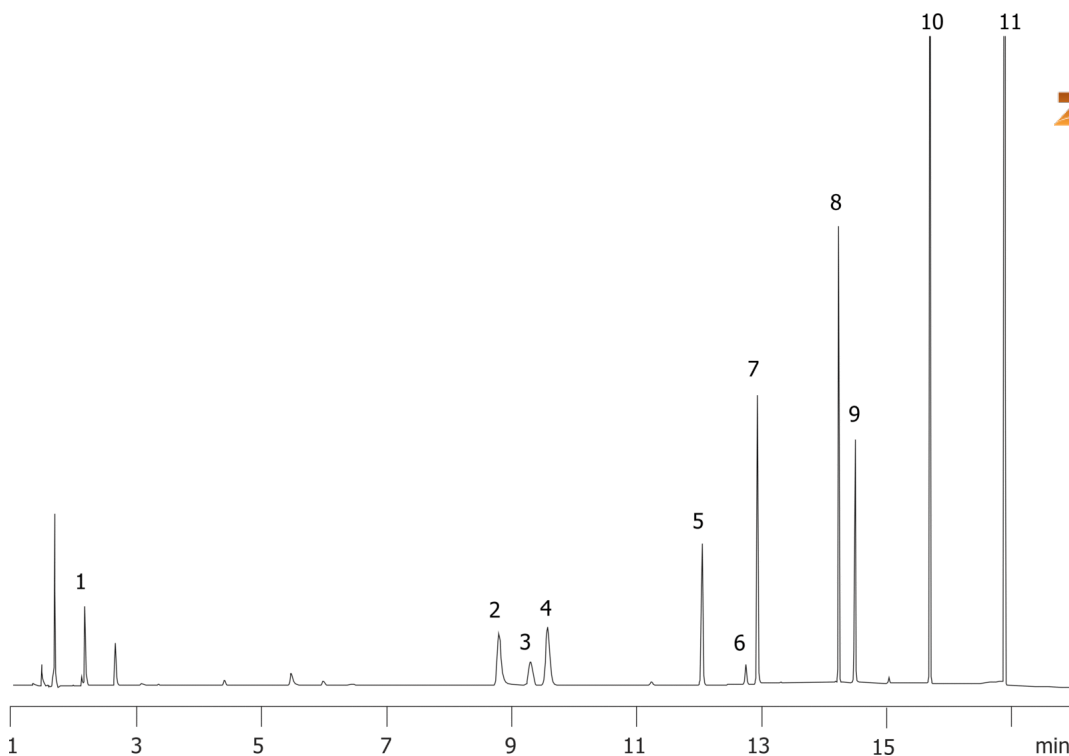


EPA Method 552.2: Haloacetic Acids on ZB-1701 30x0.25x0.25**Column:** Zebron[™] ZB-1701, GC Cap. Column 30 m x 0.25 mm x 0.25 μ m, Ea**Phase:** 14% Cyanopropylphenyl 86% Dimethylpolysiloxane**Dimensions:** 30 meters x 0.25 mm x 0.25 μ m**Order No:** 7HG-G006-11**Oven Profile:** 45°C for 10 min to 145°C at 15°C/min for 5 min.**Carrier Gas:** Constant Flow Helium, 1.9 mL/min**Injection:** split 20:1 1 μ L @ 200°C**Detection:** Electron Capture (ECD) (260°C)

14810

**Products used in this application:**

EPA Method 552.2: Haloacetic Acids on ZB-1701 30x0.25x0.25

ANALYTES:

- 1 Chloroacetic acid
- 2 Bromoacetic acid
- 3 Dalopan
- 4 Dichloroacetic acid
- 5 Trichloroacetic acid
- 6 1,2,3-Trichloropropane
- 7 Bromochloroacetic acid
- 8 Bromodichloroacetic acid
- 9 Dibromoacetic acid
- 10 Chlorodibromoacetic acid
- 11 Tribromoacetic acid

