


Amphetamines on ZB-5 30x0.25x0.25**Column:** Zebron[™] ZB-5, GC Cap. Column 30 m x 0.25 mm x 0.25 µm, Ea**Phase:** 5% Phenyl 95% Dimethylpolysiloxane**Dimensions:** 30 meters x 0.25 mm x 0.25 µm**Order No:** 7HG-G002-11**Oven Profile:** 70°C to 120°C at 10°C/min to 280°C at 20°C/min for 2 min**Carrier Gas:** Constant Flow Helium, 1.1 mL/min**Injection:** Split 10:1 1 µL @ 225°C**Detection:** Mass Selective (MSD) (180°C)

14818

**Products used in this application:****ANALYTES:**

- 1 Methamphetamine
- 2 Amphetamine (PFP)
- 3 Methamphetamine-PFP

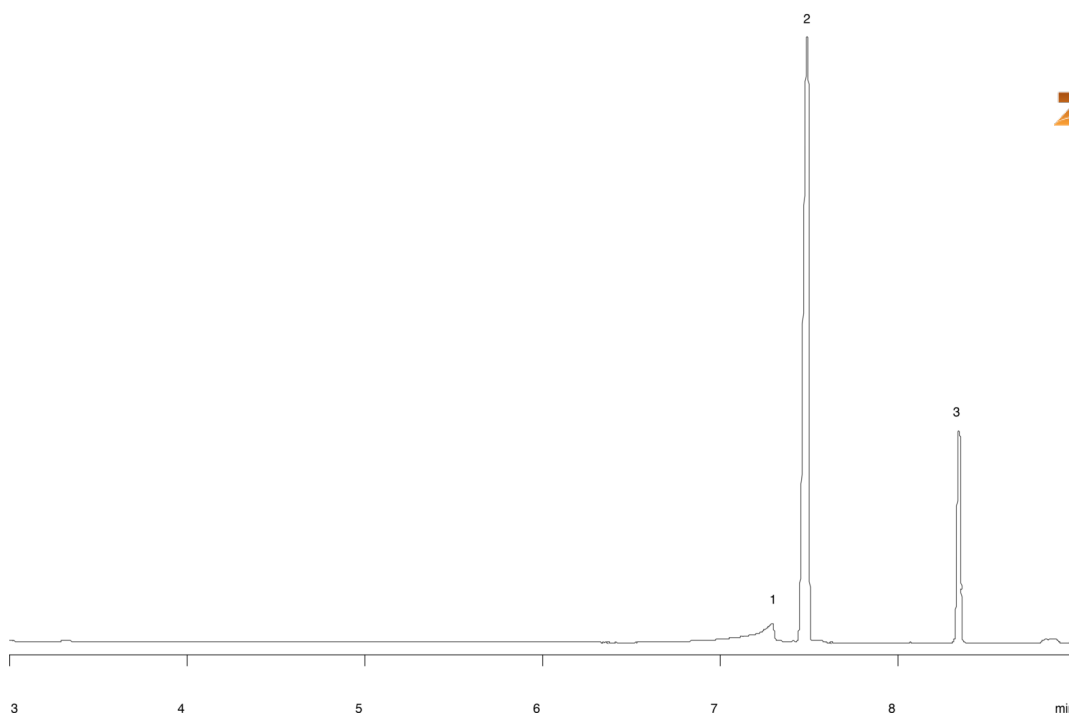