


Propane Purity using ZB-1 60x0.32x3.0**Column:** Zebron[™] ZB-1, GC Cap. Column 60 m x 0.32 mm x 3.00 µm, Ea**Phase:** 100% Dimethylpolysiloxane**Dimensions:** 60 meters x 0.32 mm x 3 µm**Order No:** 7KM-G001-36**Oven Profile:** 40°C isothermal**Carrier Gas:** Constant Flow Hydrogen, 6.5 mL/min**Injection:** split 10:1 5 µL @ 200°C**Detection:** Flame Ionization (FID) (200°C)

14820

**Products used in this application:****ANALYTES:**

- 1 Methane
- 2 Ethane
- 3 Propane
- 4 Isobutane
- 5 Butane

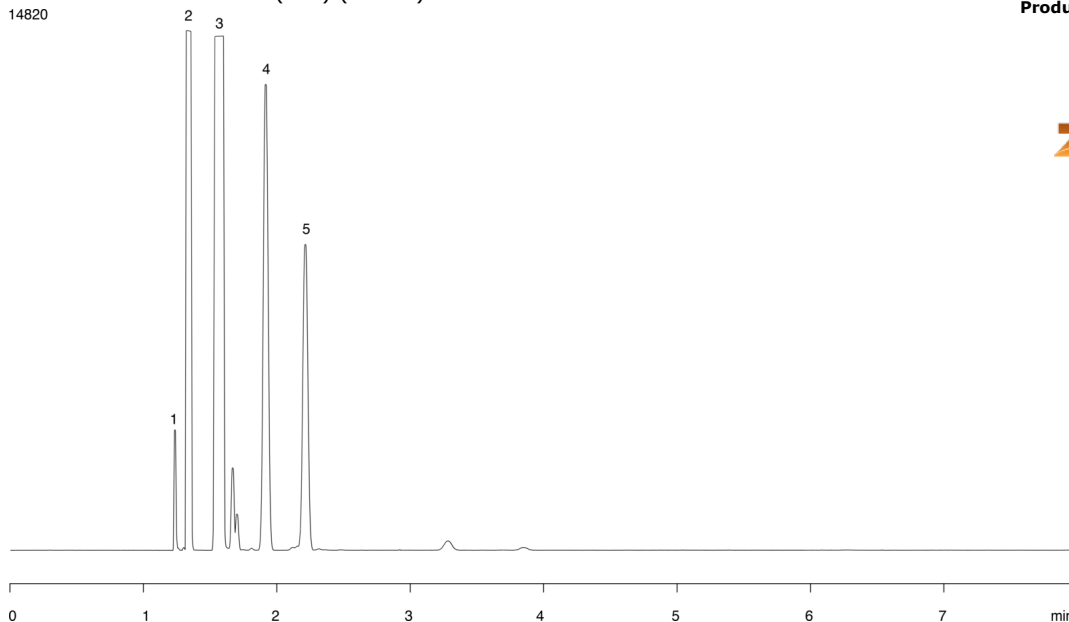