

Paraffin Wax on ZB-1

Column: Zebron™ ZB-1, GC Cap. Column 15 m x 0.53 mm x 0.15 µm, Ea

Phase: 100% Dimethylpolysiloxane

Dimensions: 15 meters x 0.53 mm x 0.15 µm

Order No: 7EK-G001-05


Oven Profile: 100°C to 360°C at 10°C/min for 20min.

Carrier Gas: Constant Flow Hydrogen, 6.7 mL/min

Injection: Splitless :1 0.2 µL @ 350°C

Detection: Flame Ionization (FID) (430°C)

14825

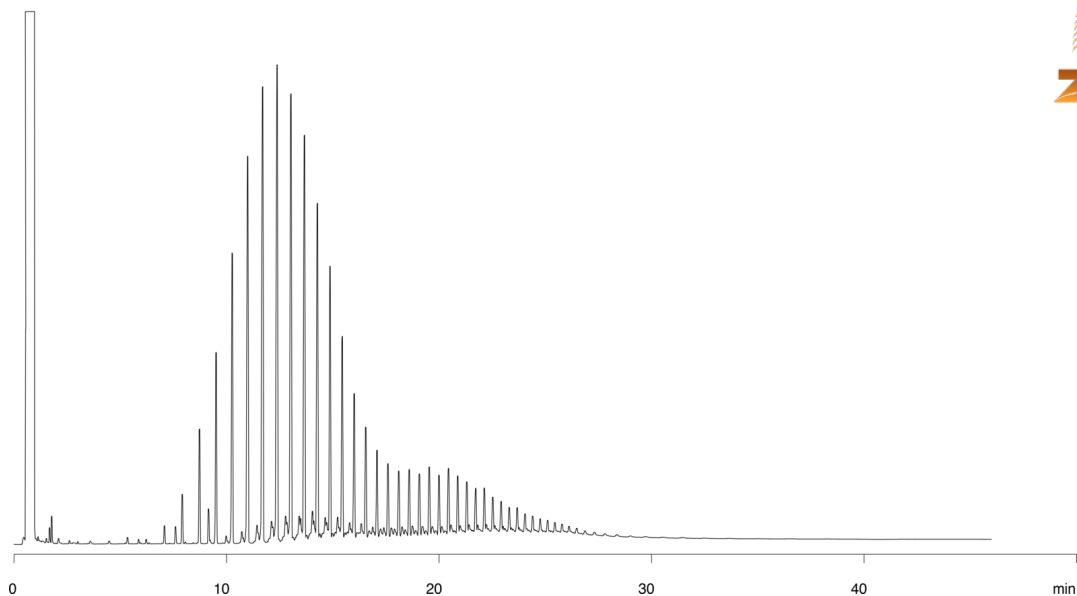

Products used in this application:



ANALYTES:

- 1 Paraffin wax

