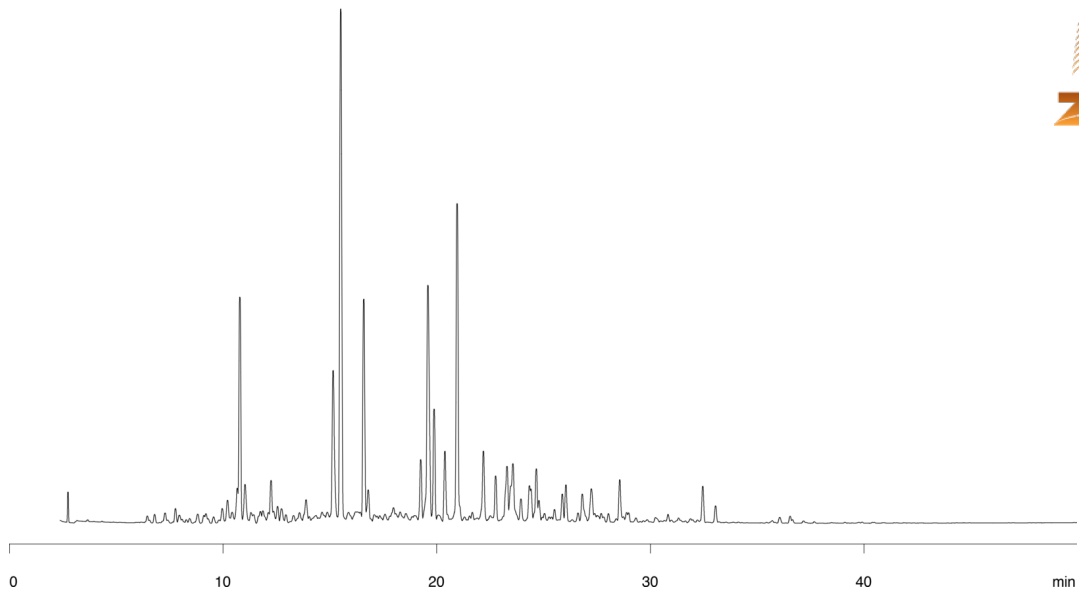


Weathered Gasoline on ZB-5

Column: Zebron™ ZB-5, GC Cap. Column 30 m x 0.53 mm x 1.50 µm, Ea
Phase: 5% Phenyl 95% Dimethylpolysiloxane
Dimensions: 30 meters x 0.53 mm x 1.5 µm
Order No: 7HK-G002-28
Oven Profile: 40°C for 5 min to 300°C at 4°C/min for 5 min.

Carrier Gas: Constant Flow Helium, 4.1 mL/min
Injection: split 30:1 1 µL @ 275°C
Detection: Flame Ionization (FID) (300°C)

14827

**ANALYTES:**

1 Gasoline

**Products used in this application:**