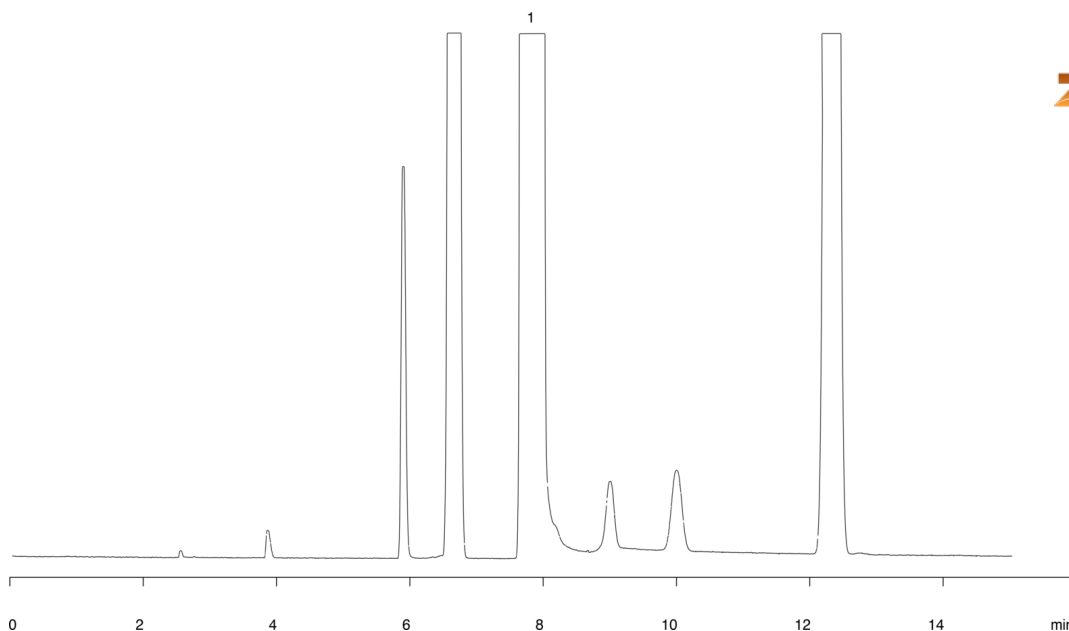


Pentane Purityl on ZB-1**Column:** Zebron[™] ZB-1, GC Cap. Column 60 m x 0.32 mm x 3.00 µm, Ea**Phase:** 100% Dimethylpolysiloxane**Dimensions:** 60 meters x 0.32 mm x 3 µm**Order No:** 7KM-G001-36**Oven Profile:** 40°C isothermal**Carrier Gas:** Constant Flow Helium, 3.5 mL/min**Injection:** split 10:1 1 µL @ 200°C**Detection:** Flame Ionization (FID) (200°C)

14830

**Products used in this application:****ANALYTES:****1** Pentane