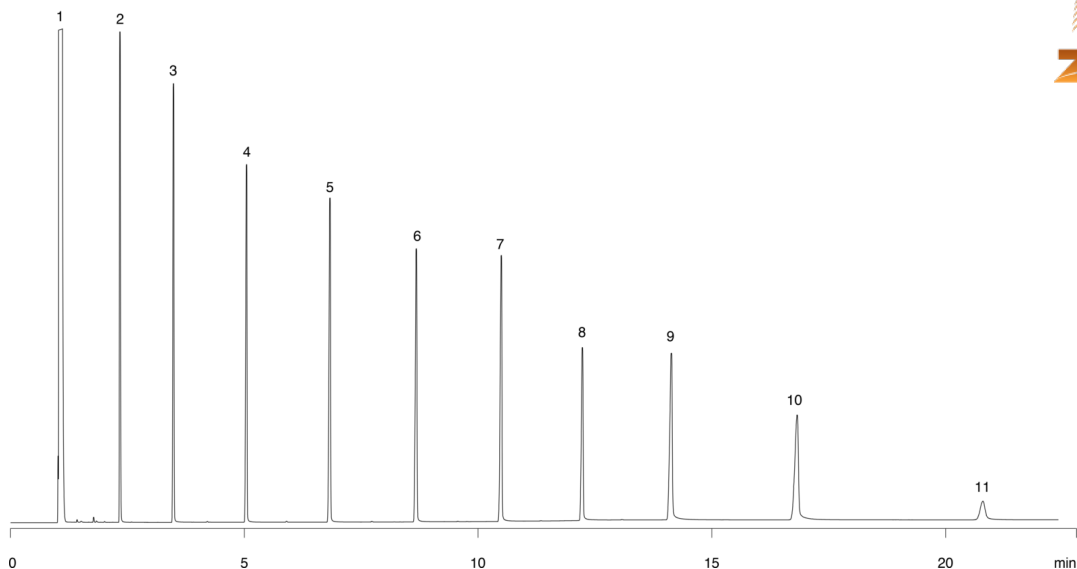


Even Carbon Free Fatty Acids on ZB-FFAP**Column:** Zebron[™] ZB-FFAP, GC Cap. Column 30 m x 0.25 mm x 0.25 μ m, Ea**Phase:** Nitroterephthalic Acid Modified Polyethylene Glycol**Dimensions:** 30 meters x 0.25 mm x 0.25 μ m**Order No:** 7HG-G009-11**Oven Profile:** 160°C to 260°C at 8°C/min**Carrier Gas:** Constant Flow Helium, 2.4 mL/min**Injection:** Split 10:1 1 μ L @ 280°C**Detection:** Flame Ionization (FID) (280°C)

14835

**Products used in this application:**

Even Carbon Free Fatty Acids on ZB-FFAP

ANALYTES:

- 1 solvent
- 2 Caproic acid
- 3 Caprylic acid
- 4 Capric acid
- 5 Lauric acid
- 6 Myristic acid
- 7 Palmitic acid
- 8 Stearic acid
- 9 Arachidic acid
- 10 Behenic acid
- 11 Linoceric acid

