

Antihistamines using Synergi Fusion-RP**Column:** SynergiTM 4 μ m Fusion-RP 80 Å, LC Column 150 x 4.6 mm, Ea**Dimensions:** 150 x 4.6 mm ID**Order No:** 00F-4424-E0**Elution Type:** Isocratic**Eluent A:** 20mm Potassium Phosphate, pH 2.5**Eluent B:** Methanol

Gradient	Step No.	Time (min)	Pct A	Pct B
Profile:	1	0	70	30

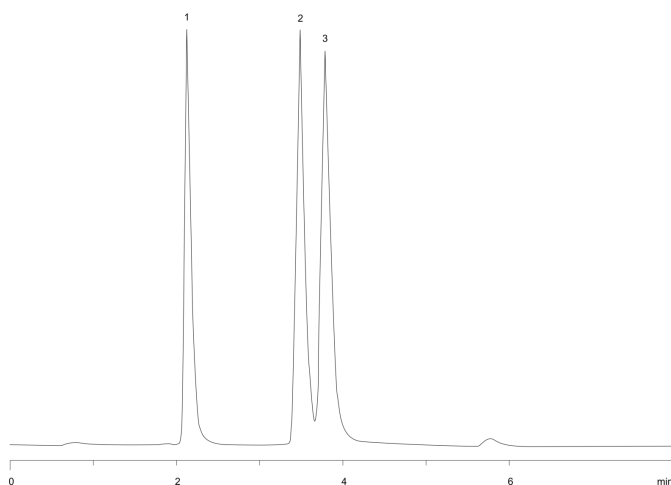
Flow Rate: 1 mL/min**Col. Temp.:** ambient**Detection:** UV-Vis Abs.-Variable Wave.(UV) @ 210 nm (ambient)**Analyst Note:** SecurityGuardTM Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, Fusion-RP 4 x 3.0mm, 10/Pk Part No.: AJ0-7557

- Holder Part No.: KJ0-4282



14839

**ANALYTES:**

- 1 Phenylephrine
- 2 Phenylpropanolamine
- 3 Pseudoephedrine

