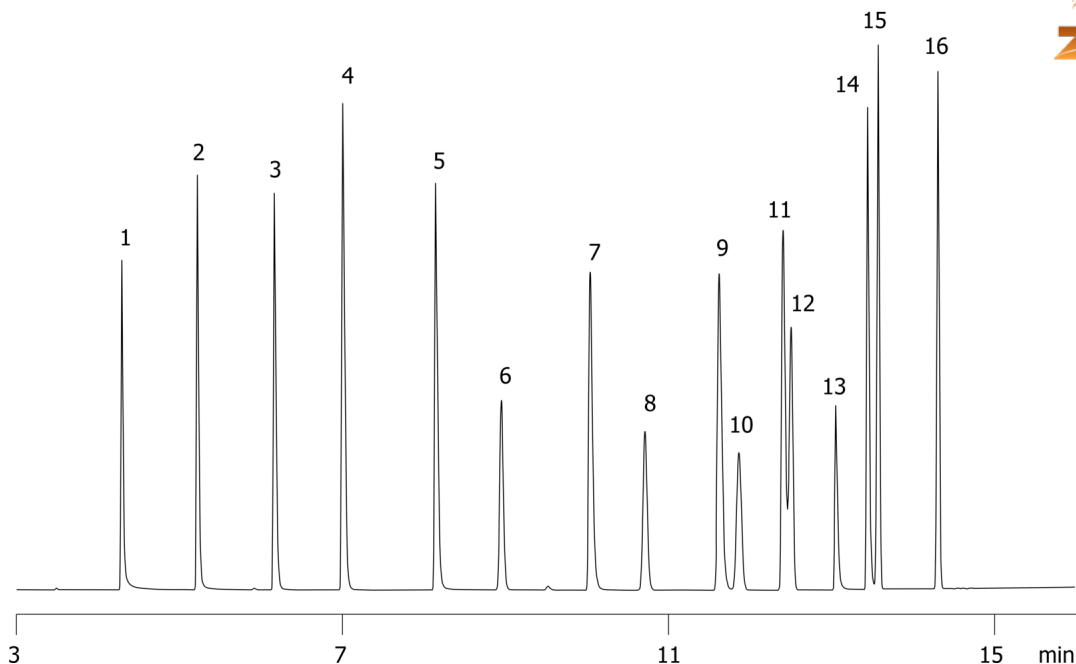


Petroleum Oxygenates ZB-1 60x0.25x1.0 ASTM 4815

Column: Zebtron™ ZB-1, GC Cap. Column 60 m x 0.25 mm x 1.00 µm, Ea
Phase: 100% Dimethylpolysiloxane
Dimensions: 60 meters x 0.25 mm x 1 µm
Order No: 7KG-G001-22
Oven Profile: 40C for 6 min to 50°C at 5°C/min for 4 min to 175°C at 25°C/min for 5 min.
Carrier Gas: Constant Flow Helium, 1.8 mL/min
Injection: Split 40:1 0.2 µL @ 200°C
Detection: Flame Ionization (FID) (250°C)
14849

**Products used in this application:**

Petroleum Oxygenates ZB-1 60x0.25x1.0 ASTM 4815

ANALYTES:

- 1 Methanol
- 2 Ethanol
- 3 Isopropanol
- 4 tert-Butanol
- 5 1-Propanol
- 6 Methyl tert-butyl ether
- 7 sec-Butanol
- 8 Isopropyl ether
- 9 Isobutanol
- 10 tert-Butyl ethyl ether
- 11 tert-Pentanol
- 12 Methylcyclopentane
- 13 1,2-Dimethoxyethane
- 14 n-Butyl alcohol
- 15 Benzene
- 16 tert-Amyl methyl ether

