


Nitrosamines ZB-WAX 30x0.32x0.5**Column:** Zebron[™] ZB-WAX, GC Cap. Column 30 m x 0.32 mm x 0.50 μ m, Ea**Phase:** Polyethylene Glycol (PEG)**Dimensions:** 30 meters x 0.32 mm x 0.5 μ m**Order No:** 7HM-G007-17**Oven Profile:** 40°C to 260°C at 20°C/min for 5 min.**Carrier Gas:** Constant Flow Helium, 1.4 mL/min**Injection:** split 50:1 1 μ L @ 200°C**Detection:** Mass Selective (MSD) (240°C)

14854

**Products used in this application:****ANALYTES:**

- 1 N-Nitrosodimethylamine
- 2 N-Nitrosomethylethylamine
- 3 N-Nitrosodiethylamine
- 4 N-Nitrosodi-n-propylamine
- 5 N-Nitrosodi-n-butylamine
- 6 N-Nitrosopiperidine
- 7 N-Nitrosopyrrolidine
- 8 N-Nitrosomorpholine
- 9 N-Nitrosodiphenylamine

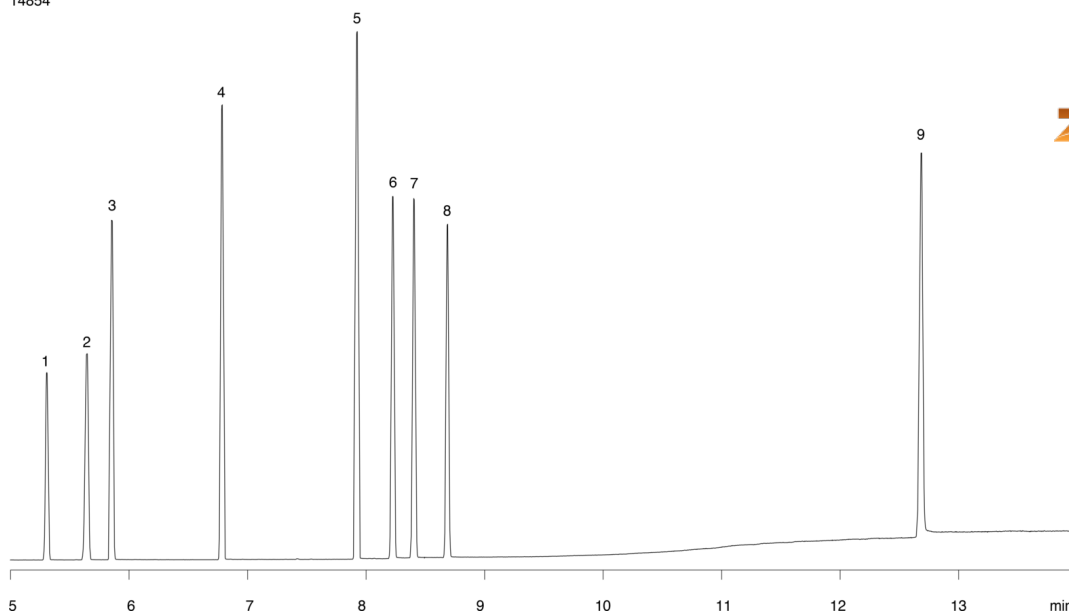