


Aniline - Phenol on ZB-5 30x0.32x1.0**Column:** Zebtron[™] ZB-5, GC Cap. Column 30 m x 0.32 mm x 1.00 µm, Ea**Phase:** 5% Phenyl 95% Dimethylpolysiloxane**Dimensions:** 30 meters x 0.32 mm x 1 µm**Order No:** 7HM-G002-22**Oven Profile:** 130°C to 160°C at 5°C/min to 240°C at 20°C/min.**Carrier Gas:** Constant Flow Helium, 1.5 mL/min**Injection:** split 15:1 0.2 µL @ 200°C**Detection:** Flame Ionization (FID) (250°C)

14862

**Products used in this application:****ANALYTES:**

- 1 Phenol
- 2 Aniline

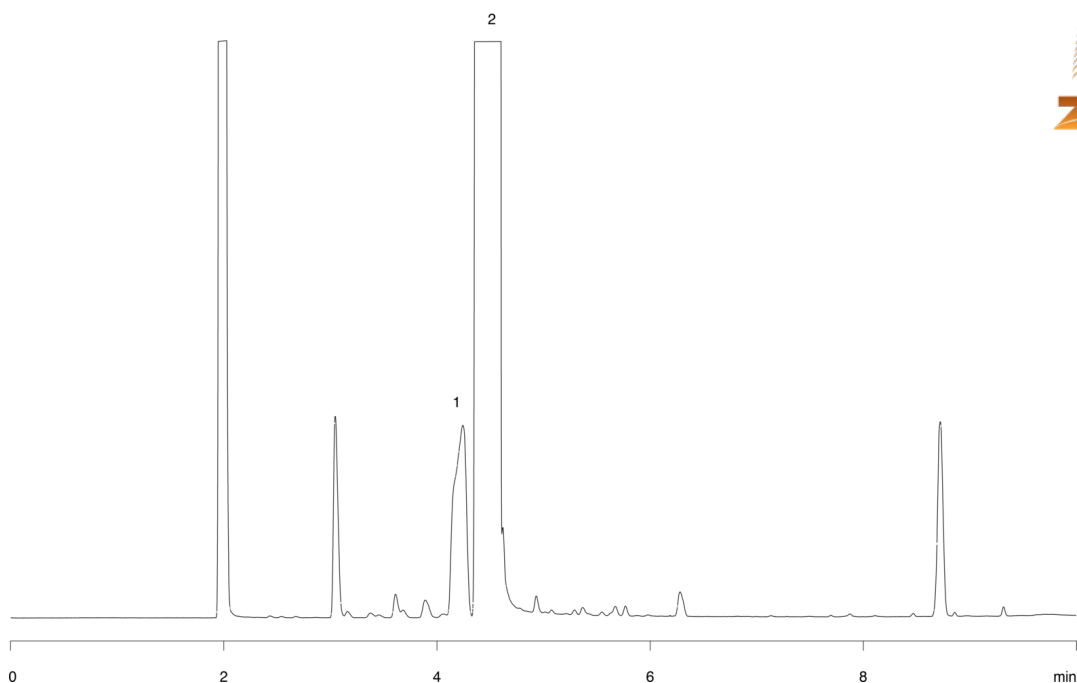