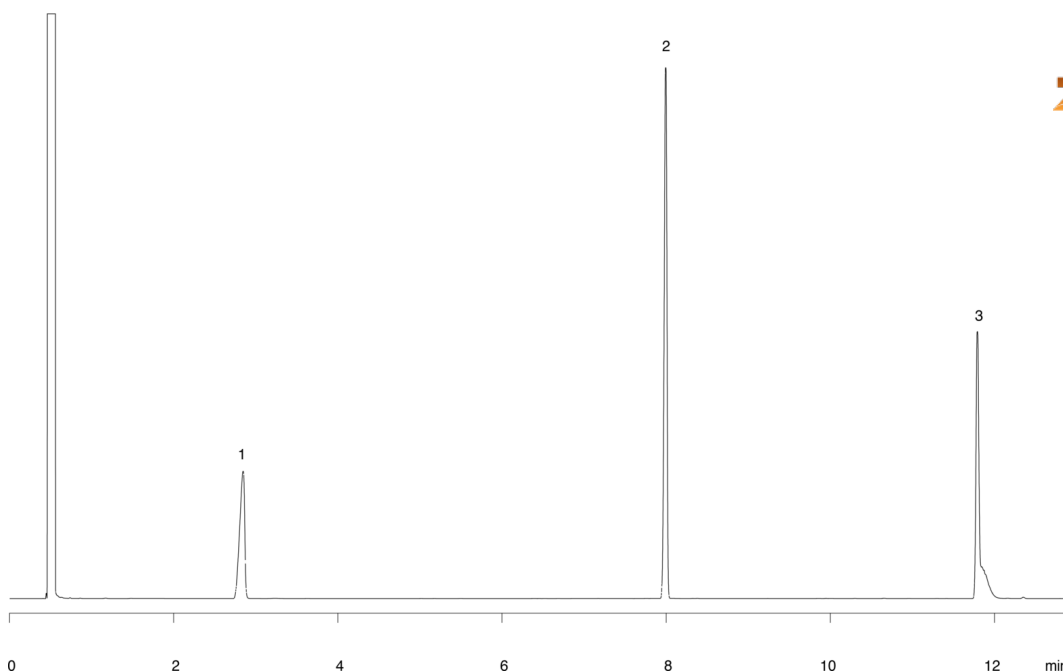


Nitrosamines - EPA Method 607 on ZB-5 15x0.53x1.5

Column: Zebtron™ ZB-5, GC Cap. Column 15 m x 0.53 mm x 1.50 µm, Ea
Phase: 5% Phenyl 95% Dimethylpolysiloxane
Dimensions: 15 meters x 0.53 mm x 1.5 µm
Order No: 7EK-G002-28
Oven Profile: 40°C for 4 min to 240°C at 20°C/min for 5 min.

Carrier Gas: Constant Flow Helium, 11.4 mL/min
Injection: split 15:1 1 µL @ 240°C
Detection: Flame Ionization (FID) (250°C)

14865

**Products used in this application:****ANALYTES:**

- 1 N-Nitrosodimethylamine
- 2 N-Nitrosodi-n-propylamine
- 3 N-Nitrosodiphenylamine

