


Nitroaromatics - EPA Method 609: ZB-5 15x0.53x1.5**Column:** Zebron™ ZB-5, GC Cap. Column 15 m x 0.53 mm x 1.50 µm, Ea**Phase:** 5% Phenyl 95% Dimethylpolysiloxane**Dimensions:** 15 meters x 0.53 mm x 1.5 µm**Order No:** 7EK-G002-28**Oven Profile:** 100°C to 330°C at 15°C/min for 5 min.**Carrier Gas:** Constant Flow Helium, 3 mL/min**Injection:** split 15:1 1 µL @ 250°C**Detection:** Flame Ionization (FID) (300°C)

14866

**Products used in this application:****ANALYTES:**

- 1 Nitrobenzene
- 2 Isophorone
- 3 2,6-Dinitrotoluene
- 4 2,4-Dinitrotoluene

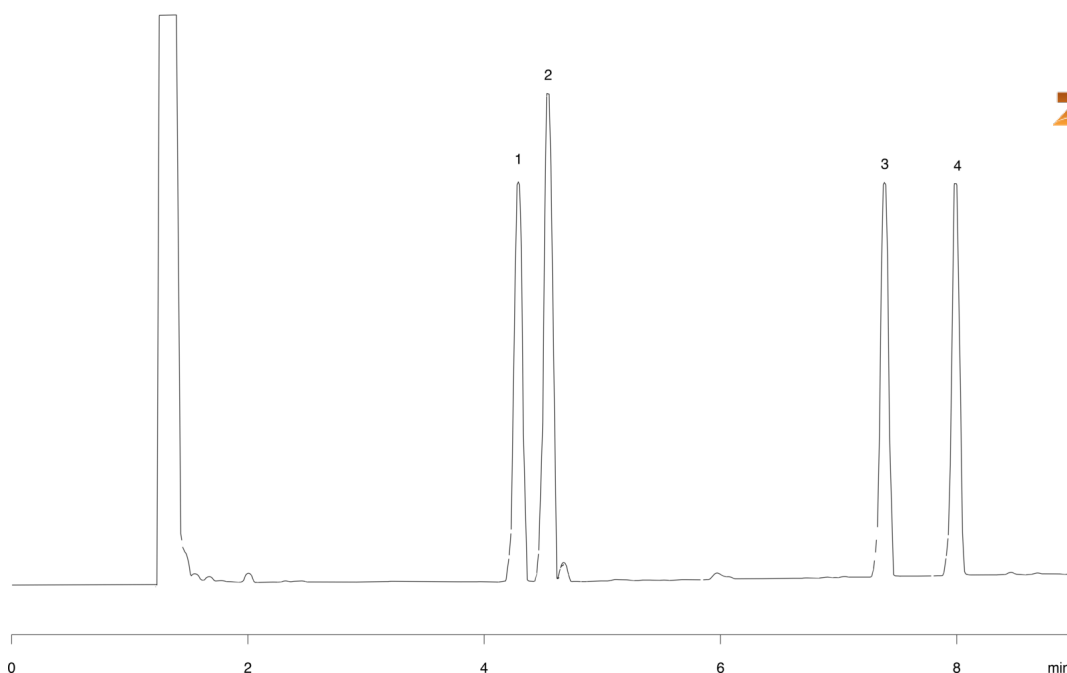