


Chlorinated Hydrocarbons: EPA Method 612 - ZB-5 30x0.53x0.5**Column:** Zebron™ ZB-5, GC Cap. Column 30 m x 0.53 mm x 0.50 µm, Ea**Phase:** 5% Phenyl 95% Dimethylpolysiloxane**Dimensions:** 30 meters x 0.53 mm x 0.5 µm**Order No:** 7HK-G002-17**Oven Profile:** 40°C to 250°C at 8°/min for 5 min.**Carrier Gas:** Constant Flow Helium, 4.9 mL/min**Injection:** split 15:1 1 µL @ 250°C**Detection:** Electron Capture (ECD) (300°C)

14879

**Products used in this application:****ANALYTES:**

- 1 1,3-Dichlorobenzene
- 2 1,4-Dichlorobenzene
- 3 1,2-Dichlorobenzene
- 4 Hexachloroethane
- 5 1,2,4-Trichlorobenzene
- 6 Hexachlorobutadiene
- 7 Hecachlorocyclopentadiene
- 8 2-Chloronaphthalene
- 9 Hexachlorobenzene

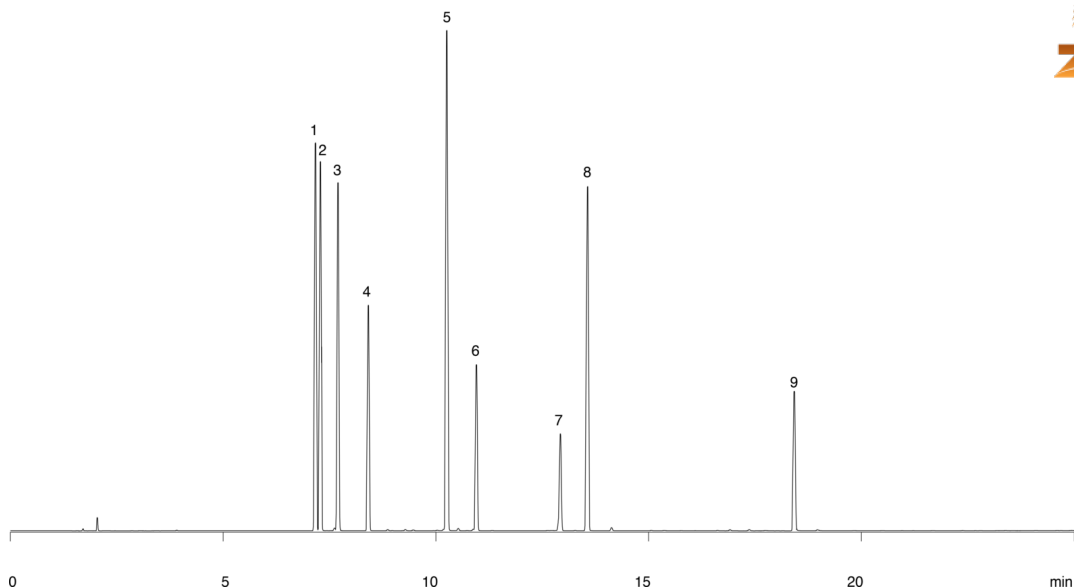