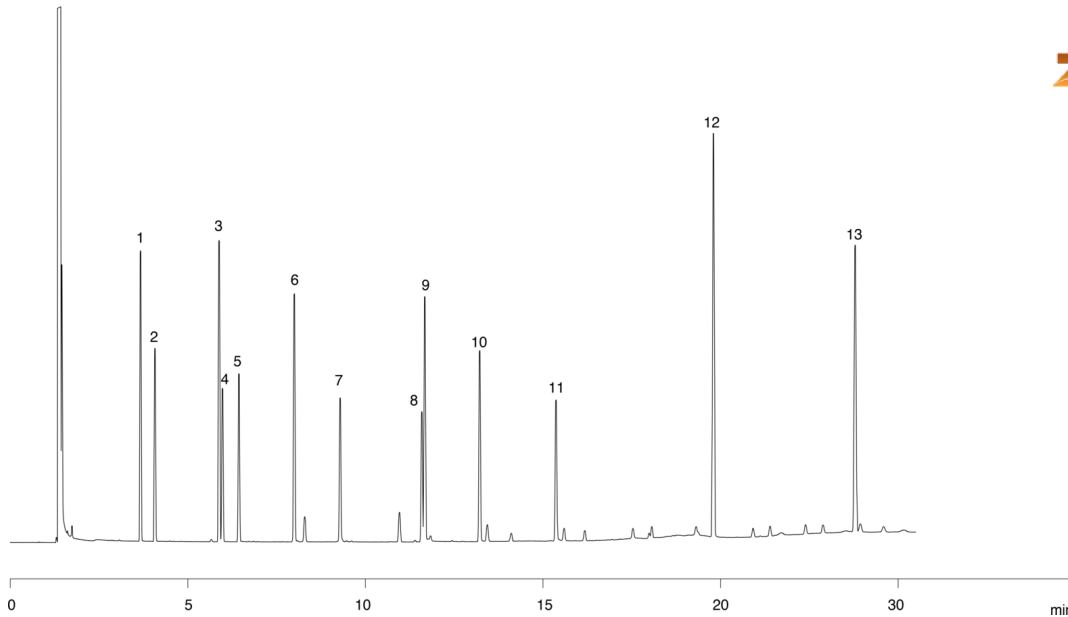


Phenols and Benzidines: EPA Method 604 + 605 - ZB-5 30x0.32x1.0**Column:** Zebron™ ZB-5, GC Cap. Column 30 m x 0.32 mm x 1.00 µm, Ea**Phase:** 5% Phenyl 95% Dimethylpolysiloxane**Dimensions:** 30 meters x 0.32 mm x 1 µm**Order No:** 7HM-G002-22**Oven Profile:** 110°C to 290°C at 8°C/min for 3 min**Carrier Gas:** Constant Flow Helium, 2.2 mL/min**Injection:** split 20:1 1 µL @ 290°C**Detection:** Flame Ionization (FID) (310°C)

14886

**Products used in this application:**

Phenols and Benzidines: EPA Method 604 + 605 - ZB-5 30x0.32x1.0

ANALYTES:

- 1 Phenol
- 2 2-Chlorophenol
- 3 2-Nitrophenol
- 4 2,4-Dimethylphenol
- 5 2,4-Dichlorophenol
- 6 4-Chloro-3-methylphenol
- 7 2,4,6-Trichlorophenol
- 8 2,4-Dinitrophenol
- 9 4-Nitrophenol
- 10 2-Methyl-4,6-dinitrophenol
- 11 Pentachlorophenol
- 12 Benzidine
- 13 3,3'-Dichlorobenzidine

