

Aqueous Stability using Synergi Fusion-RP

Column: SynergiTM 4 μ m Fusion-RP 80 Å, LC Column 150 x 4.6 mm, Ea

Dimensions: 150 x 4.6 mm ID

Order No: 00F-4424-E0

Elution Type: Isocratic

Eluent A: 20mM Potassium Phosphate buffer, pH 2.5

Gradient	Step No.	Time (min)	Pct A
Profile:	1	0	100

Flow Rate: 1 mL/min

Col. Temp.: 30 °C

Detection: UV-Vis Abs.-Variable Wave.(UV) @ 254 nm (ambient)



Analyst Note: SecurityGuardTM Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, Fusion-RP 4 x 3.0mm, 10/Pk Part No.: AJ0-7557

- Holder Part No.: KJ0-4282



14891

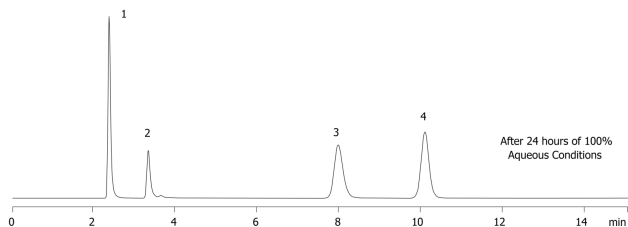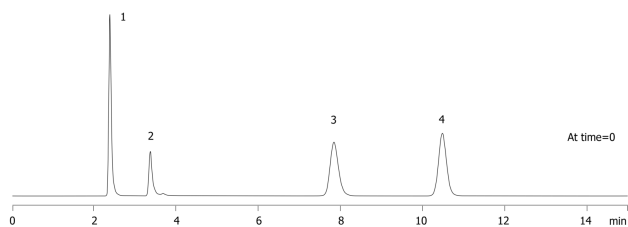

ANALYTES:

- 1 Thiourea
- 2 Adenine
- 3 Guanosine-5-monophosphate
- 4 Thymine

