


Odd Carbon Free Fatty Acids - ZBFFAP 30x0.25x0.25**Column:** Zebron[™] ZB-FFAP, GC Cap. Column 30 m x 0.25 mm x 0.25 μ m, Ea**Phase:** Nitroterephthalic Acid Modified Polyethylene Glycol**Dimensions:** 30 meters x 0.25 mm x 0.25 μ m**Order No:** 7HG-G009-11**Oven Profile:** 160°C to 260°C at 8°C/min for 15 min**Carrier Gas:** Constant Flow Helium, 2.4 mL/min**Injection:** Split 40:1 0.2 μ L @ 280°C**Detection:** Flame Ionization (FID) (280°C)

14896

**Products used in this application:****ANALYTES:**

- 1** solvent
- 2** Nonanoic acid
- 3** Undecanoic acid
- 4** Tridecanoic acid
- 5** Pentadecanoic acid
- 6** Heptadecanoic acid
- 7** Nonadecanoic acid
- 8** Heneicosanoic acid
- 9** Tricosanoic acid

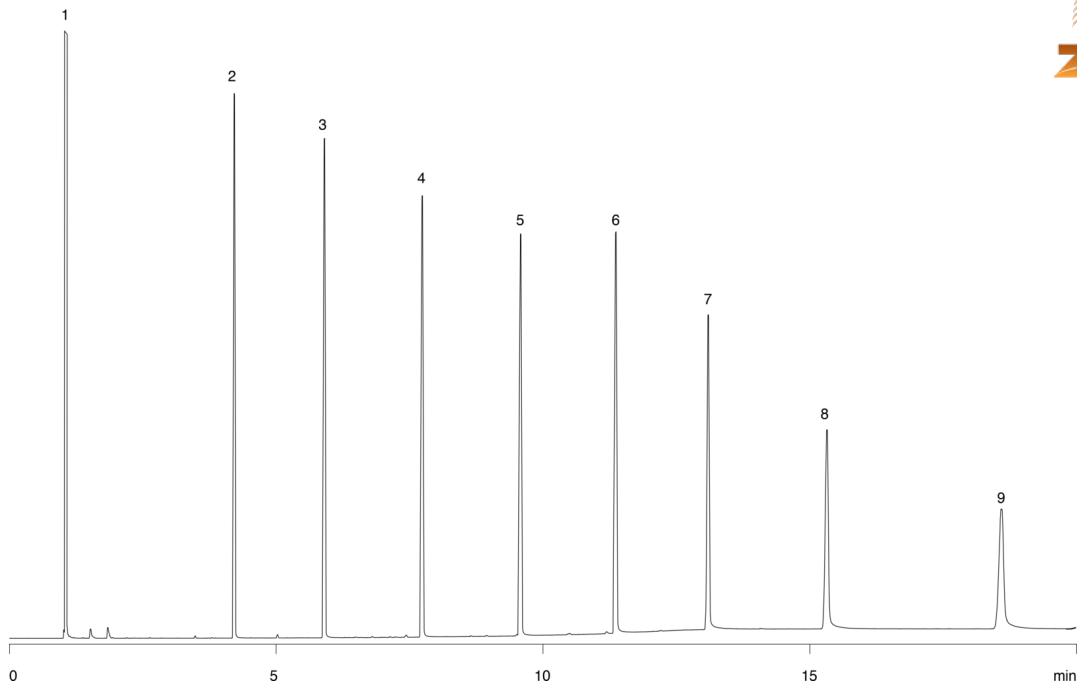