


Flavors Analysis ZB-WAX 30x0.25x0.25

Column: Zebron™ ZB-WAX, GC Cap. Column 30 m x 0.25 mm x 0.25 µm, Ea**Phase:** Polyethylene Glycol (PEG)**Dimensions:** 30 meters x 0.25 mm x 0.25 µm**Order No:** 7HG-G007-11**Oven Profile:** 50°C to 250°C at 6°C/min for 3 min.**Carrier Gas:** Constant Flow Helium, 1 mL/min**Injection:** Split 100:1 1 µL @ 250°C**Detection:** Mass Selective (MSD) (275°C)**Analyst Note:** Test Mix Courtesy of Flavour Chemistry Laboratory Frutarom (UK) Limited

Products used in this application:

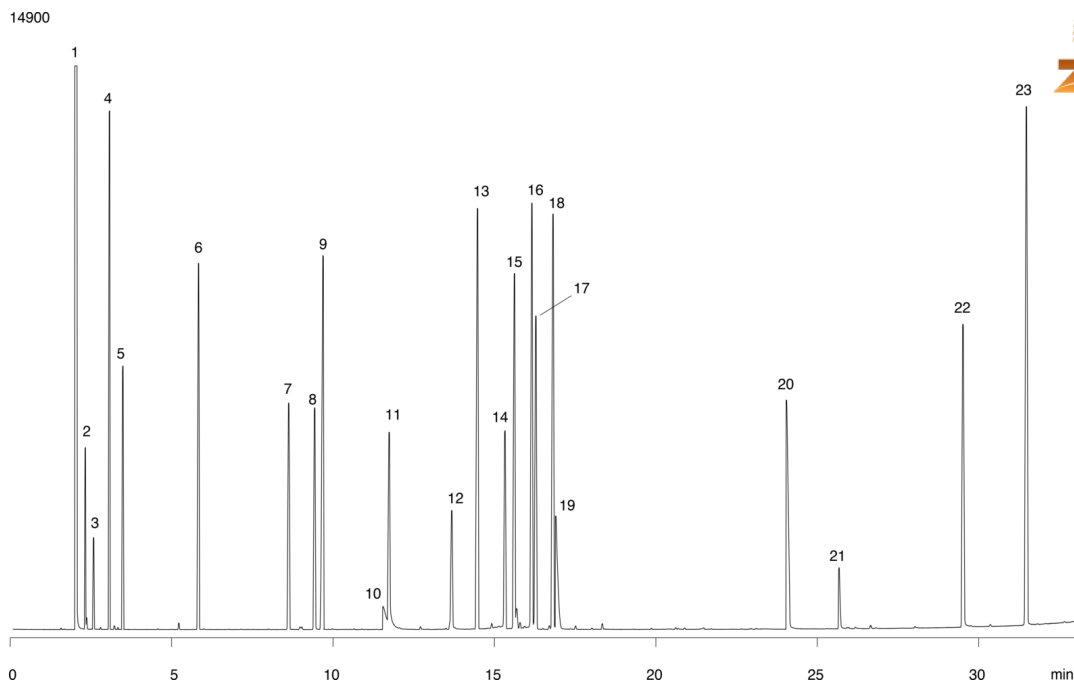

Flavors Analysis ZB-WAX 30x0.25x0.25

ANALYTES:

- 1 Acetone
- 2 Ethyl acetate
- 3 Ethanol
- 4 Decane
- 5 Ethyl butyrate
- 6 Limonene
- 7 2,3-Dimethylpyrazine
- 8 (Z)-3-Hexenol
- 9 Tetradecane
- 10 Acetic acid
- 11 Decanal
- 12 Propylene glycol
- 13 Ethyl decanoate
- 14 Neral
- 15 alpha-Terpineol
- 16 Neryl Acetate
- 17 Geranial
- 18 Decanol
- 19 Valeric acid
- 20 Nonanoic acid
- 21 Decanoic acid
- 22 Vanillin
- 23 Anthracene

