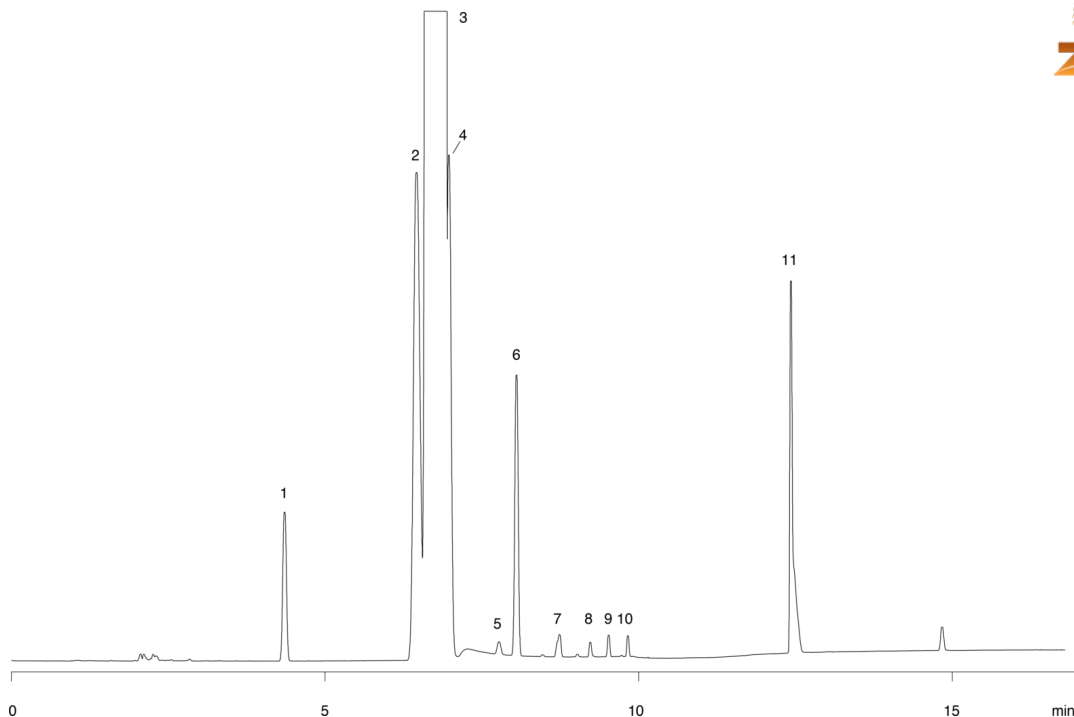


p-Xylene Purity on ZB-WAX**Column:** Zebron[™] ZB-WAX, GC Cap. Column 30 m x 0.53 mm x 1.00 μ m, Ea**Phase:** Polyethylene Glycol (PEG)**Dimensions:** 30 meters x 0.53 mm x 1 μ m**Order No:** 7HK-G007-22**Oven Profile:** 70°C for 7 min to 190 at 25°C/min for 5 min.**Carrier Gas:** Constant Flow Helium, 4 mL/min**Injection:** split 30:1 1 μ L @ 225°C**Detection:** Flame Ionization (FID) (325°C)

14941

**Products used in this application:**

p-Xylene Purity on ZB-WAX

ANALYTES:

- 1 Toluene
- 2 Ethylbenzene
- 3 p-Xylene
- 4 m-Xylene
- 5 Cumene
- 6 o-Xylene
- 7 1-Ethyl-2-methylbenzene
- 8 Styrene
- 9 Trimethylbenzene
- 10 Diethylbenzene
- 11 Methylbenzaldehyde

