


Acid Test Mix on ZB-5MS**Column:** Zebron[™] ZB-5MS, GC Cap. Column 30 m x 0.25 mm x 0.25 μ m, Ea**Phase:** 5% Phenyl-Arylene 95% Dimethylpolysiloxane**Dimensions:** 30 meters x 0.25 mm x 0.25 μ m**Order No:** 7HG-G010-11**Oven Profile:** 50°C for 2 min to 250°C at 15°C/min for 5 min**Carrier Gas:** Constant Flow Helium, 1.15 mL/min**Injection:** Split :1 2 μ L @ 250°C**Detection:** Flame Ionization (FID) (325°C)

14951

**Products used in this application:****ANALYTES:**

- 1 2,3-Butanediol
- 2 2-Ethylhexanoic acid
- 3 2,4-Dimethylphenol
- 4 Benzoic acid
- 5 Tridecane
- 6 Tetradecane
- 7 Pentadecane

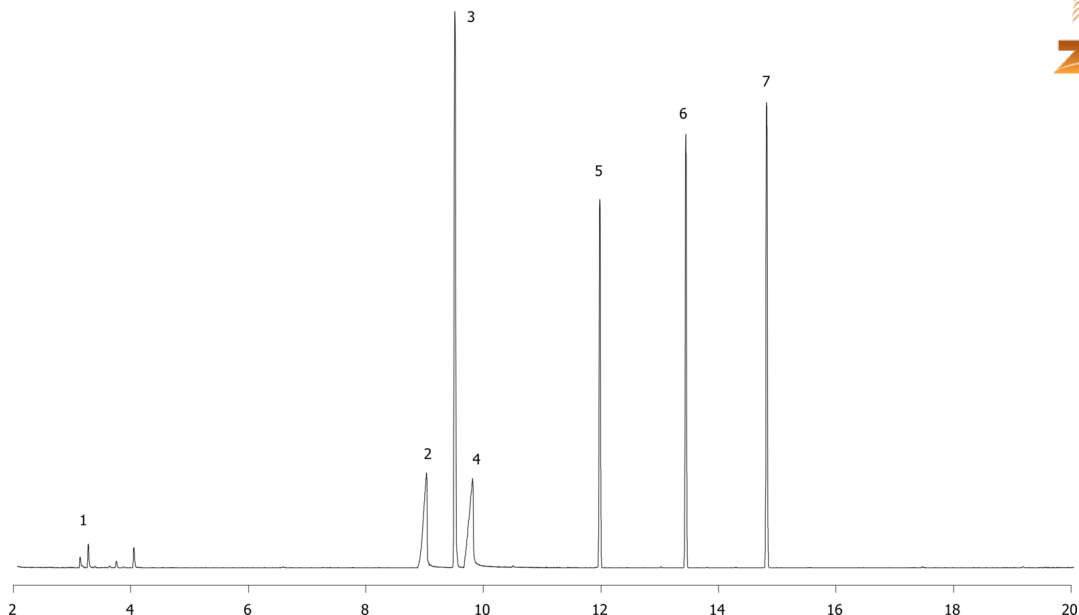