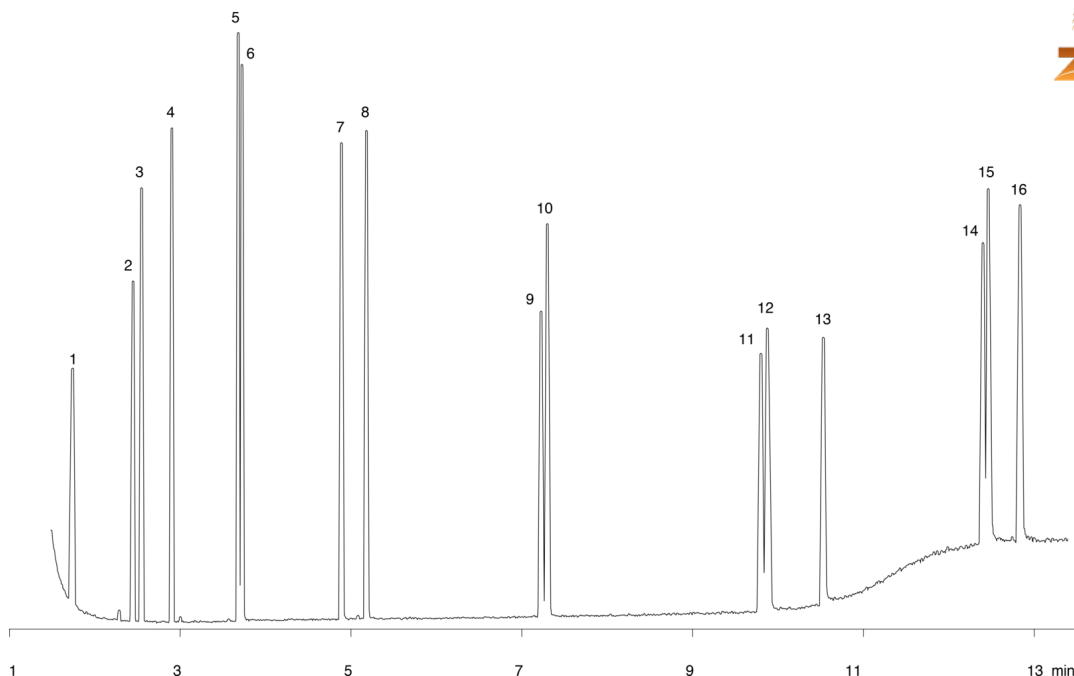


PAHs on ZB-5MS in under 13 minutes

Column: Zebtron[™] ZB-5MS, GC Cap. Column 30 m x 0.25 mm x 0.25 µm, Ea
Phase: 5% Phenyl-Arylene 95% Dimethylpolysiloxane
Dimensions: 30 meters x 0.25 mm x 0.25 µm
Order No: 7HG-G010-11
Oven Profile: 155°C to 250°C at 25°C/min to 280°C at 5°C/min to 320°C at 25°C/min for 2 min
Carrier Gas: Constant Flow Helium, 1.4 mL/min
Injection: split 15:1 1 µL @ 285°C
Detection: Mass Selective (MSD) (240°C)

14952

**Products used in this application:**

PAHs on ZB-5MS in under 13 minutes

ANALYTES:

- 1 Naphthalene
- 2 Acenaphthylene
- 3 Acenaphthene
- 4 Fluorene
- 5 Phenanthrene
- 6 Anthracene
- 7 Fluoranthene
- 8 Pyrene
- 9 Benzo[a]anthracene
- 10 Chrysene
- 11 Benzo[b]fluoranthene
- 12 Benzo[k]fluoranthene
- 13 Benzo[a]pyrene
- 14 Indeno[1,2,3-cd]pyrene
- 15 Dibenzo[a,h]anthracene
- 16 Benzo[g,h,i]perylene

