

## Oligosaccharides using Rezex RSO

**Column:** Rezex™ RSO-Oligosaccharide Ag+ 4%, LC Column 200 x 10 mm, Ea

**Dimensions:** 200 x 10 mm ID

**Order No:** 00P-0133-N0

**Elution Type:** Isocratic

**Eluent A:** Water

| Gradient | Step No. | Time (min) | Pct A |
|----------|----------|------------|-------|
| Profile: | 1        | 0          | 100   |

**Flow Rate:** 0.3 mL/min

**Col. Temp.:** 75 °C

**Detection:** Evaporative Light Scattering @ nm (ambient)


14964

**REZEX™**  
HPLC COLUMNS for Carbohydrates, Oligosaccharides, and Organic Acid Analysis



Products used in this application:

**REZEX™**

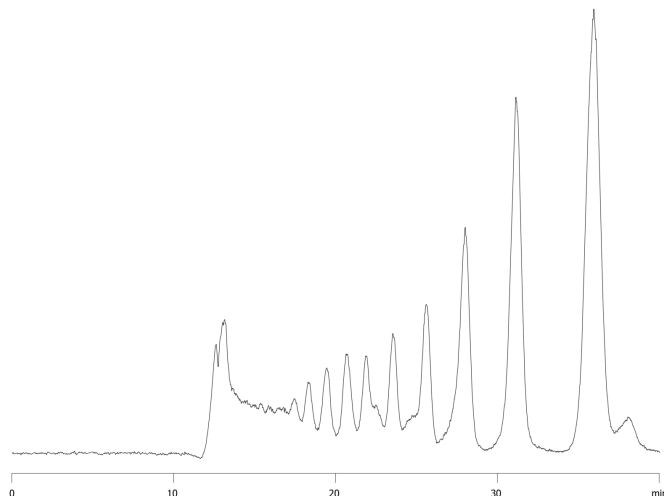

### ANALYTES:

1 Malto-oligosaccharides

