

Saccharides using Rezex RAM**Column:** N/A**Dimensions:** N/A**Order No:** N/A**Elution Type:** N/A**Eluent A:** Water

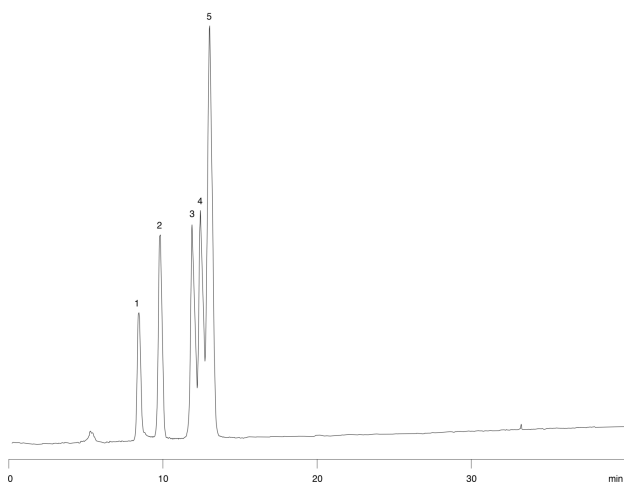
Gradient	Step No.	Time (min)	Pct A
Profile:	1	0	100

Flow Rate: 0.6 mL/min**Col. Temp.:** 75 °C**Detection:** Evaporative Light Scattering @ nm (ambient)**Analyst Note:** SecurityGuardTM Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, Carbo-Ag 4 x 3.0mm, 10/Pk Part No.: AJ0-4491

- Holder Part No.: KJ0-4282

14966

**ANALYTES:**

- 1 Melezitose
- 2 Mannose
- 3 Maltose
- 4 Fructose
- 5 Glucose

