

Wood saccharides using Rezex RPM

Column: Rezex™ RPM-Monosaccharide Pb+2 (8%), LC Column 300 x 7.8 mm, Ea

Dimensions: 300 x 7.8 mm ID

Order No: 00H-0135-K0

Elution Type: Isocratic

Eluent A: Water

Gradient	Step No.	Time (min)	Pct A
Profile:	1	0	100

Flow Rate: 0.6 mL/min

Col. Temp.: 85 °C

Detection: Evaporative Light Scattering @ nm (ambient)

Analyst Note: SecurityGuard™ Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, Carbo-Pb 4 x 3.0mm, 10/Pk Part No.: AJ0-4492

- Holder Part No.: KJ0-4282


14967

REZEX™
HPLC COLUMNS for Carbohydrates, Oligosaccharides, and Organic Acid Analysis



Products used in this application:

REZEX™

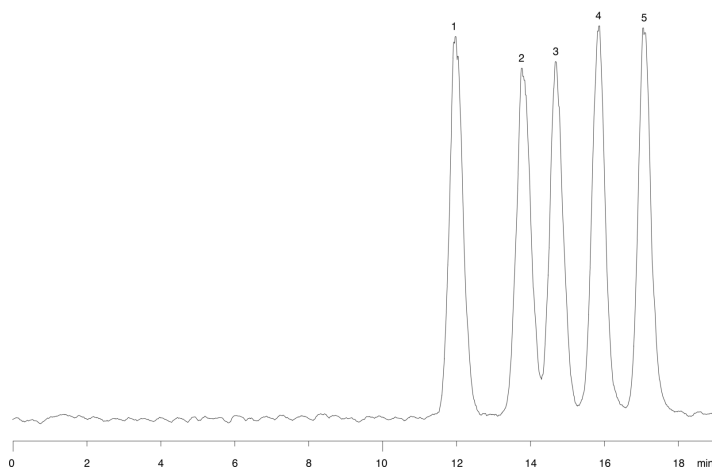

ANALYTES:

- 1 D-Cellobiose
- 2 alpha-D-Glucose
- 3 D-(+)-Xylose
- 4 D-(+)-Galactose
- 5 D-(-)-Arabinose

