

Carbohydrates analysis using Rezex RAM

Column: N/A

Dimensions: N/A

Order No: N/A

Elution Type: N/A

Eluent A: Water

Gradient	Step No.	Time (min)	Pct A
Profile:	1	0	100

Flow Rate: 0.6 mL/min

Col. Temp.: 80 °C


Detection: Evaporative Light Scattering @ nm (ambient)

Analyst Note: SecurityGuard™ Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, Carbo-Ag 4 x 3.0mm, 10/Pk Part No.: AJ0-4491

- Holder Part No.: KJ0-4282

14974

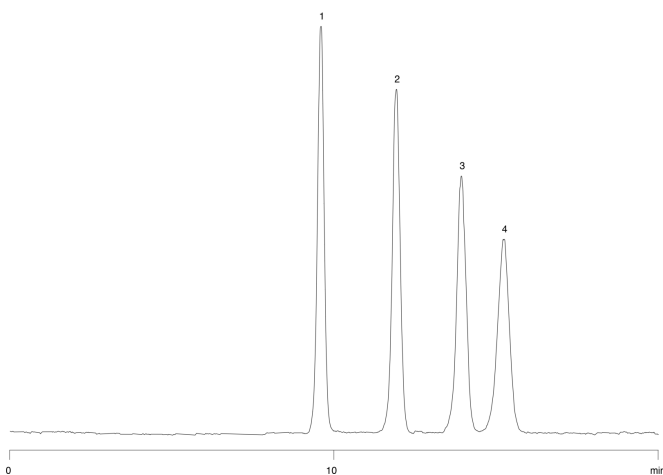

ANALYTES:

- 1 Sucrose
- 2 Glucose
- 3 Arabinose
- 4 Ribose

