

Saccharide analysis using Rezex RCM

Column: Rezex™ RCM-Monosaccharide Ca+2 (8%), LC Column 300 x 7.8 mm, Ea

Dimensions: 300 x 7.8 mm ID

Order No: 00H-0130-K0

Elution Type: Isocratic

Eluent A: Water

Gradient Profile:	Step No.	Time (min)	Pct A
1	1	0	100

Flow Rate: 0.6 mL/min

Col. Temp.: 80 °C


Detection: Evaporative Light Scattering @ nm (ambient)

Analyst Note: SecurityGuard™ Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, Carbo-Ca 4 x 3.0mm, 10/Pk Part No.: AJ0-4493

- Holder Part No.: KJ0-4282

14975

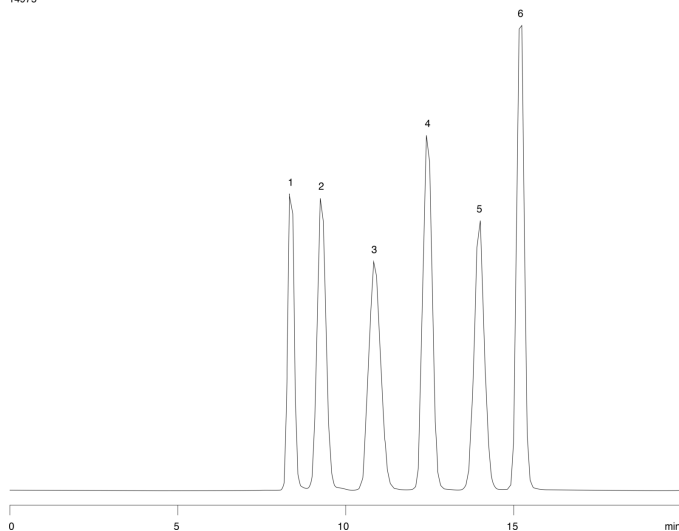

ANALYTES:

- 1 Melezitose
- 2 Maltose
- 3 Glucose
- 4 Mannose
- 5 Fructose
- 6 Ribitol

REZEX™
HPLC COLUMNS for Carbohydrates, Oligosaccharides, and Organic Acid Analysis



Products used in this application:

REZEX™

