

Polar Acids

Column: Gemini® 5 µm C18 110 Å, LC Column 150 x 4.6 mm, Ea

Dimensions: 150 x 4.6 mm ID

Order No: 00F-4435-E0

Elution Type: Isocratic

Eluent A: 20mM Potassium Phosphast, pH 6.0

Eluent B: Acetonitrile

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	70	30

Flow Rate: 1 mL/min

Col. Temp.: 25 °C

Detection: UV-Vis Abs.-Variable Wave.(UV) @ 202 nm (25 °C)

Analyst Note: SecurityGuard™ Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, Gemini C18 4 x 3.0mm, 10/Pk Part No.: AJ0-7597


- Holder Part No.: KJ0-4282



Products used in this application:



15000

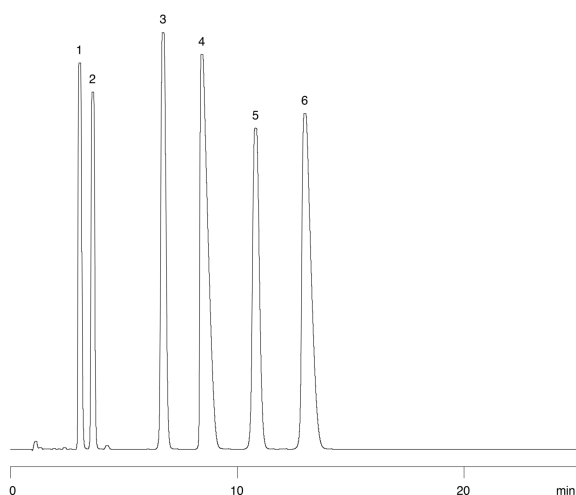

ANALYTES:

- 1** Tolmetin
Retention Time: 3.046 min
- 2** Naproxen
Retention Time: 3.627 min
- 3** Fenoprofen
Retention Time: 6.728 min
- 4** Diflunisal
Retention Time: 8.443 min
- 5** Indomethacin
Retention Time: 10.808 min
- 6** Ibuprofen
Retention Time: 12.988 min

