

Caffeine & Metabolites, pH 4.5

Column: Gemini® 5 µm C18 110 Å, LC Column 150 x 4.6 mm, Ea

Dimensions: 150 x 4.6 mm ID

Order No: 00F-4435-E0

Elution Type: Gradient

Eluent A: 20mM Ammonium acetate, pH 4.5

Eluent B: Acetonitrile

Gradient	Step No.	Time (min)	Pct A	Pct B
Profile:	1	0	98	2
	2	15	80	20

Flow Rate: 1 mL/min

Col. Temp.: 25 °C

Detection: UV-Vis Abs.-Variable Wave.(UV) @ 273 nm (25 °C)

Analyst Note: SecurityGuard™ Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, Gemini C18 4 x 3.0mm, 10/Pk Part No.: AJ0-7597


- Holder Part No.: KJ0-4282



Products used in this application:



15013

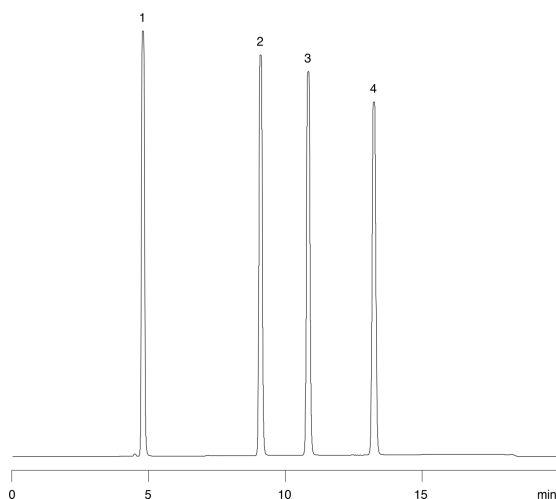

ANALYTES:

1 Xanthine

Retention Time: 4.791 min

2 Theobromine

Retention Time: 9.092 min

3 Theophylline

Retention Time: 10.833 min

4 Caffeine

Retention Time: 13.235 min

