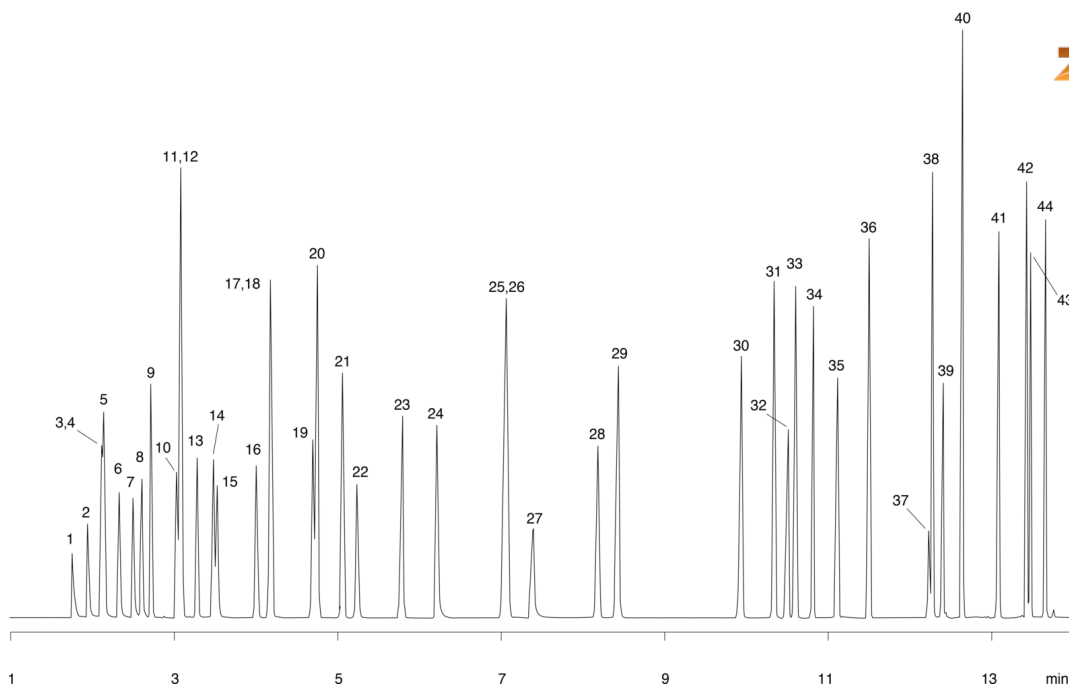


Multiple Solvents Mix 2 on ZB-5MS**Column:** Zebron[™] ZB-5MS, GC Cap. Column 30 m x 0.25 mm x 0.50 µm, Ea**Phase:** 5% Phenyl-Arylene 95% Dimethylpolysiloxane**Dimensions:** 30 meters x 0.25 mm x 0.5 µm**Order No:** 7HG-G010-17**Oven Profile:** 40°C for 2 to 70°C at 4°C/min to 280°C at 25°C/min**Carrier Gas:** Constant Flow Helium, 0.9 mL/min**Injection:** split 100:1 0.2 µL @ 200°C**Detection:** Mass Selective (MSD) (240°C)

15420

**Products used in this application:**

Multiple Solvents Mix 2 on ZB-5MS**ANALYTES:**

- | | |
|-----------------------------------|--------------------------|
| 1 Methanol | 41 Benzyl alcohol |
| 2 Ethanol | 42 Acetophenone |
| 3 Acetonitrile | 43 o-Toluidine |
| 4 Acetone | 44 Nitrobenzene |
| 5 Isopropanol | |
| 6 Acrylonitrile | |
| 7 Allyl alcohol | |
| 8 Propanol | |
| 9 Methyl tert-butyl ether | |
| 10 2-Butanone | |
| 11 Hexane | |
| 12 Isopropyl ether | |
| 13 Ethyl acetate | |
| 14 Chloroform | |
| 15 Tetrahydrofuran | |
| 16 1,2-Dichloroethane | |
| 17 Benzene | |
| 18 Butanol | |
| 19 Propyl ether | |
| 20 2,2,4-Trimethylpentane | |
| 21 Trichloroethylene | |
| 22 1,4-Dioxane | |
| 23 Methyl cyclohexane | |
| 24 Pyridine | |
| 25 Pentanol | |
| 26 Toluene | |
| 27 Dimethyl formamide | |
| 28 Octane | |
| 29 Tetrachloroethylene | |
| 30 Chlorobenzene | |
| 31 Ethylbenzene | |
| 32 Hexanol | |
| 33 p-Xylene | |
| 34 n-Butyl ether | |
| 35 Cyclohexanone | |
| 36 Anisole | |
| 37 Benzaldehyde | |
| 38 Heptanol | |
| 39 Aniline | |
| 40 1,2,4-Trimethyl benzene | |

