


Butane Purity ZB1 60x0.32x0.25**Column:** ZebtronTM ZB-1, GC Cap. Column 60 m x 0.32 mm x 0.25 µm, Ea**Phase:** 100% Dimethylpolysiloxane**Dimensions:** 60 meters x 0.32 mm x 0.25 µm**Order No:** 7KM-G001-11**Oven Profile:** 40°C Isothermal**Carrier Gas:** Constant Flow Hydrogen, 6.5 mL/min**Injection:** split 10:1 5 µL @ 200°C**Detection:** Flame Ionization (FID) (200°C)

15426

**Products used in this application:****ANALYTES:**

- 1 Ethane
- 2 Propane
- 3 Isobutane
- 4 Butane

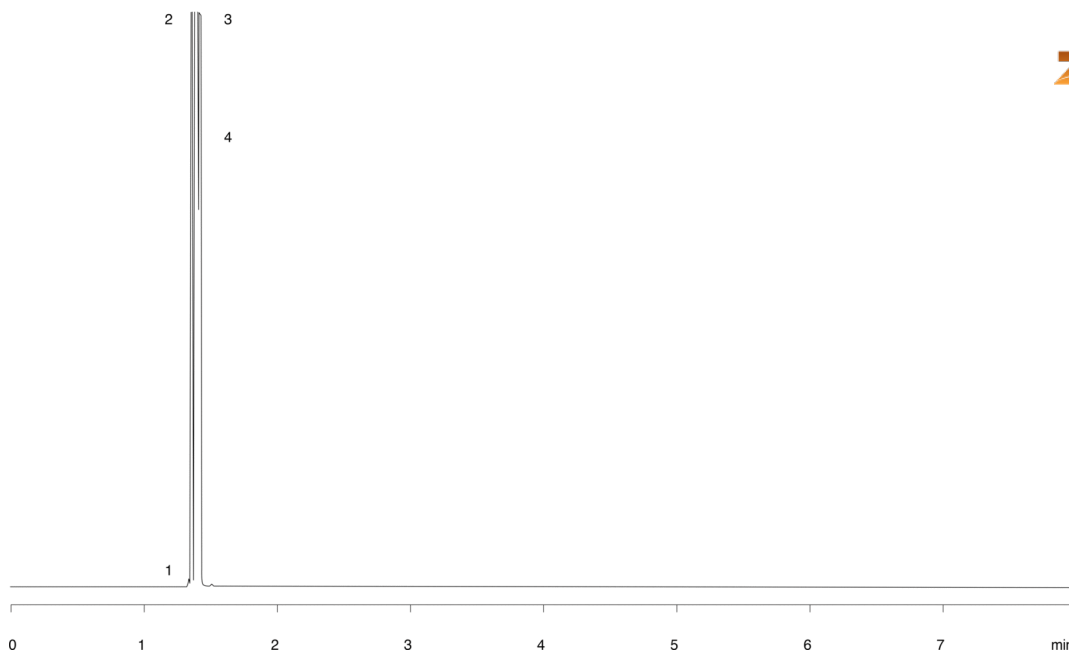