


Chlorinated Hydrocarbons on ZB-5MS**Column:** Zebron[™] ZB-5MS, GC Cap. Column 30 m x 0.25 mm x 0.25 µm, Ea**Phase:** 5% Phenyl-Arylene 95% Dimethylpolysiloxane**Dimensions:** 30 meters x 0.25 mm x 0.25 µm**Order No:** 7HG-G010-11**Oven Profile:** 60°C to 120 at 10°C/min to 230°C at 20°C/min for 2 min**Carrier Gas:** Constant Flow Helium, 1 mL/min**Injection:** split 20:1 1 µL @ 250°C**Detection:** Electron Capture (ECD) (320°C)

15428

**Products used in this application:****ANALYTES:**

- 1** 1,3-Dichlorobenzene
- 2** 1,4-Dichlorobenzene
- 3** 1,2-Dichlorobenzene
- 4** Hexachloroethane
- 5** 1,2,4-Trichlorobenzene
- 6** Hexachlorobutadiene
- 7** Hexachlorocyclopentadiene
- 8** 2-Chloronaphthalene
- 9** Hexachlorobenzene

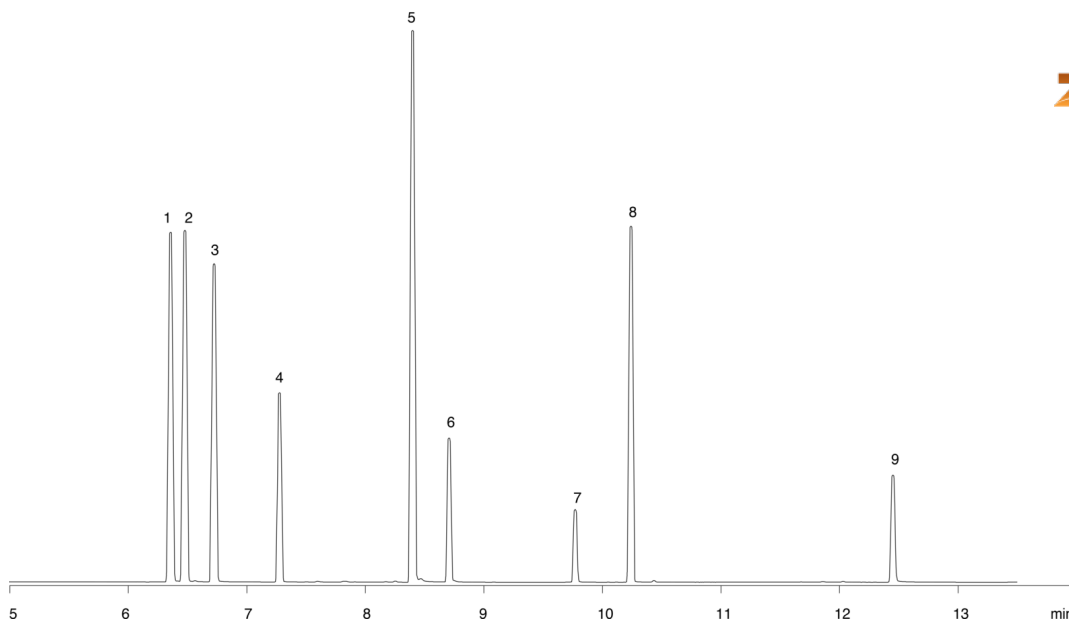