


Nitroaromatics on ZB-5MS**Column:** Zebron[™] ZB-5MS, GC Cap. Column 30 m x 0.25 mm x 0.25 µm, Ea**Phase:** 5% Phenyl-Arylene 95% Dimethylpolysiloxane**Dimensions:** 30 meters x 0.25 mm x 0.25 µm**Order No:** 7HG-G010-11**Oven Profile:** 100°C to 320°C at 15°C/min for 5 min**Carrier Gas:** Constant Flow Helium, 1 mL/min**Injection:** split 20:1 1 µL @ 250°C**Detection:** Mass Selective (MSD) (240°C)

15430

**Products used in this application:****ANALYTES:**

- 1 Nitrobenzene
- 2 Isophorone
- 3 2,6-Dinitrotoluene
- 4 2,4-Dinitrotoluene

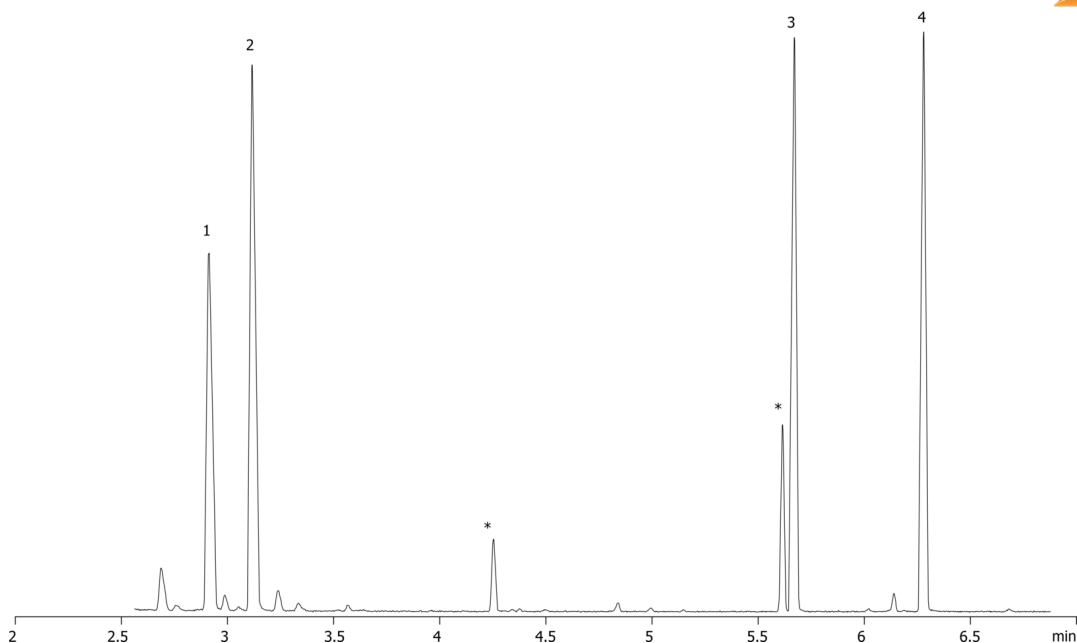