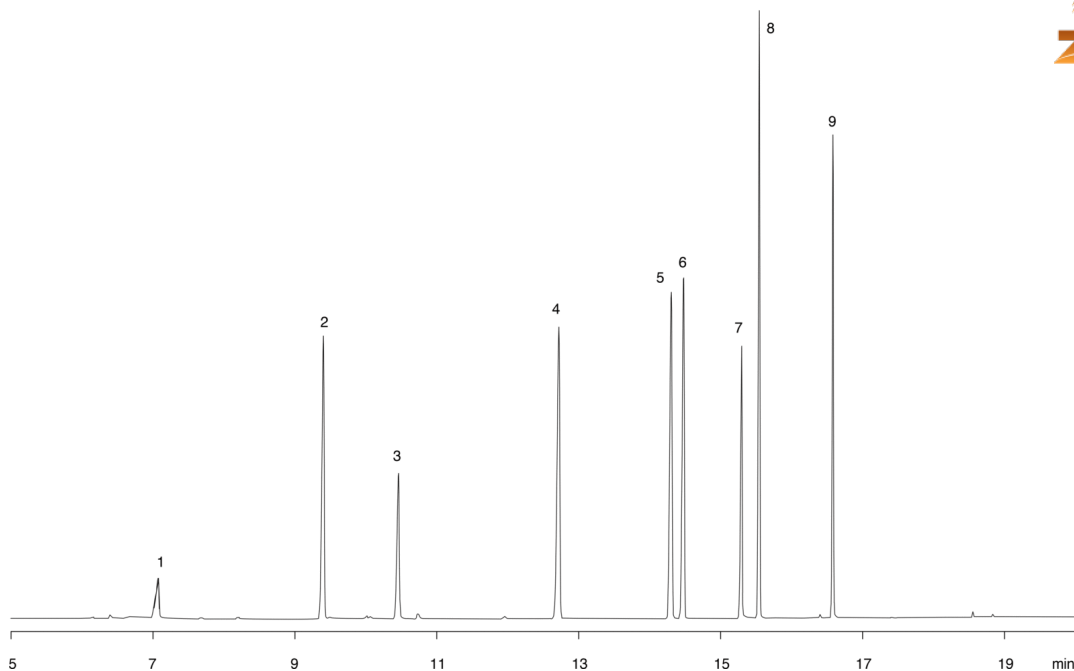


Organophosphorous Pesticides on ZB-5MS

Column: Zebron[™] ZB-5MS, GC Cap. Column 30 m x 0.25 mm x 0.25 µm, Ea
Phase: 5% Phenyl-Arylene 95% Dimethylpolysiloxane
Dimensions: 30 meters x 0.25 mm x 0.25 µm
Order No: 7HG-G010-11
Oven Profile: 100°C to 140°C at 20°C/min to 200°C at 5°C/min to 300°C at 25°C/min for 5 min
Carrier Gas: Constant Flow Helium, 1 mL/min
Injection: split 20:1 1 µL @ 250°C
Detection: Mass Selective (MSD) (240°C)

15436

**Products used in this application:****ANALYTES:**

- 1** Trichlorfon
- 2** Thionazin
- 3** Dicrotophos
- 4** Fonofos
- 5** Dichlofenthion
- 6** Chlorpyrifos-methyl
- 7** Fenitrothion
- 8** Aspon
- 9** Crotoxyphos

