

Benzodiazepines Standards on Gemini C18**Column:** Gemini® 5 µm C18 110 Å, LC Column 150 x 3 mm, Ea**Dimensions:** 150 x 3 mm ID**Order No:** 00F-4435-Y0**Elution Type:** Gradient**Eluent A:** 0.1% Formic acid in Water**Eluent B:** 0.1% Formic acid in Acetonitrile

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	80	20
	2	8	20	80
	3	9	20	80

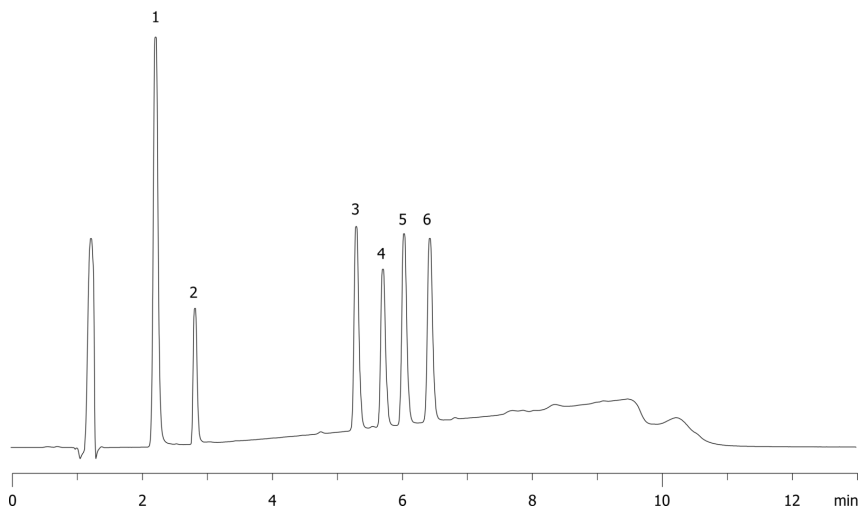
Flow Rate: 0.8 mL/min**Col. Temp.:** ambient**Detection:** UV-Vis Abs.-Variable Wave.(UV) @ 254 nm (ambient)**Analyst Note:** SecurityGuard™ Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, Gemini C18 4 x 2.0mm, 10/Pk Part No.: AJ0-7596

- Holder Part No.: KJ0-4282

**Products used in this application:**

15456

**ANALYTES:**

- 1** Chlordiazepoxide
Retention Time: 2.22 min
- 2** Flurazepam
Retention Time: 2.83 min
- 3** Oxazepam
Retention Time: 5.31 min
- 4** Clonazepam
Retention Time: 5.72 min
- 5** Temazepam
Retention Time: 6.05 min
- 6** Diazepam
Retention Time: 6.45 min

