

Beta Blockers on Onyx C18

Column: Onyx™ Monolithic C18, LC Column 100 x 4.6 mm, Ea

Dimensions: 100 x 4.6 mm ID

Order No: CH0-7643

Elution Type: Isocratic

Eluent A: 20mM Sodium Phosphate, pH 3.0

Eluent B: Methanol

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	70	30

Flow Rate: 2 mL/min


Col. Temp.: ambient

Detection: UV-Vis Abs.-Variable Wave.(UV) @ 254 nm (ambient)

15487



Products used in this application:

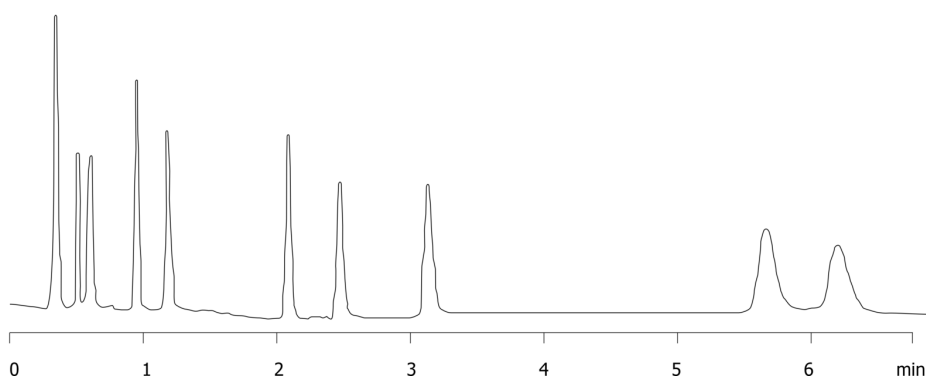

Beta Blockers on Onyx C18

ANALYTES:

- 1** Fumaric acid
- 2** Pindolol
- 3** Nadolol
- 4** Pafenolol
- 5** Metoprolol
- 6** Celiprolol
- 7** Carazolol
- 8** Bisoprolol
- 9** Alprenolol
- 10** Propranolol

