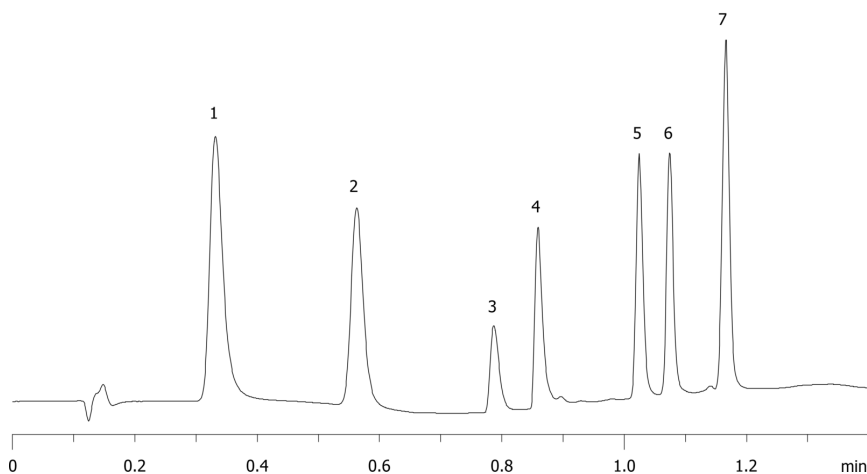


Pharmaceutical Compounds on Synergi Fusion-RP 2.5µm 20 x 2.0mm MercuryMS Cartridge**Column:** N/A**Dimensions:** N/A**Order No:** N/A**Elution Type:** N/A**Eluent A:** 0.1% Formic Acid in Water**Eluent B:** 0.1% Formic acid in Acetonitrile

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	92	8
	2	1	5	95
	3	1.01	92	8
	4	2	92	8

Flow Rate: 0.6 mL/min**Col. Temp.:** ambient**Detection:** UV-Vis Abs.-Variable Wave.(UV) @ 254 nm (ambient)

15516

**ANALYTES:**

- 1** Acetaminophen
Retention Time: 0.332 min
- 2** Caffeine
Retention Time: 0.563 min
- 3** Propanolol
Retention Time: 0.787 min
- 4** Imipramine
Retention Time: 0.859 min
- 5** Prednisolone
Retention Time: 1.024 min
- 6** Naproxen
Retention Time: 1.074 min
- 7** Indomethacine
Retention Time: 1.165 min

