

Aflatoxins using Synergi Max-RP**Column:** Synergi[™] 4 μ m Max-RP 80 Å, LC Column 150 x 3 mm, Ea**Dimensions:** 150 x 3 mm ID**Order No:** 00F-4337-Y0**Elution Type:** Isocratic**Eluent A:** 10 mM Ammonium Acetate, pH 4.5**Eluent B:** Acetonitrile

Gradient	Step No.	Time (min)	Pct A	Pct B
Profile:	1	0	70	30

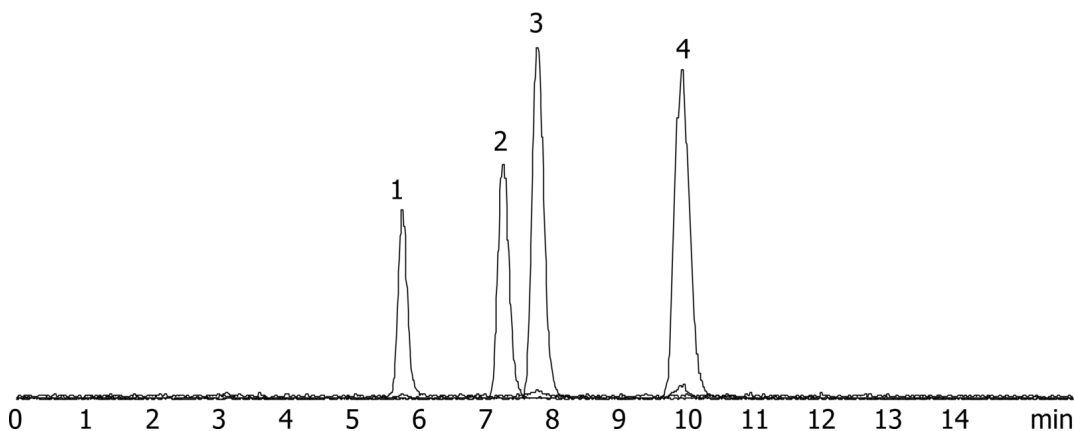
Flow Rate: 0.5 mL/min**Col. Temp.:** ambient**Detection:** Mass Spectrometer (MS) @ amu (ambient)**Analyst Note:** SecurityGuard[™] Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, C12 4 x 2.0mm, 10/Pk Part No.: AJ0-6073

- Holder Part No.: KJ0-4282



15523

**ANALYTES:**

- 1 Aflatoxin G2
- 2 Aflatoxin B2
- 3 Aflatoxin G1
- 4 Aflatoxin B1

