

Zearalone analysis using Synergi Polar-RP**Column:** SynergiTM 4 μ m Polar-RP 80 Å, LC Column 150 x 3 mm, Ea**Dimensions:** 150 x 3 mm ID**Order No:** 00F-4336-Y0**Elution Type:** Isocratic**Eluent A:** 0.1% Formic acid in Water**Eluent B:** Acetonitrile

Gradient	Step No.	Time (min)	Pct A	Pct B
Profile:	1	0	50	50

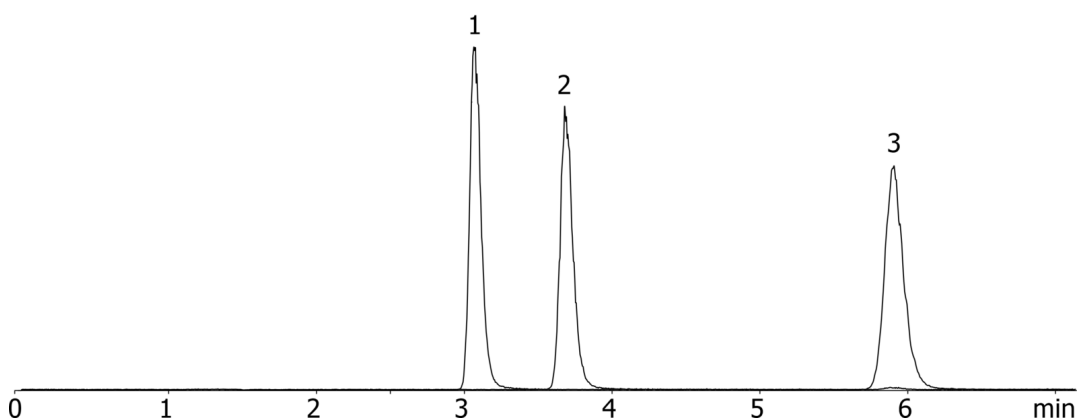
Flow Rate: 0.5 mL/min**Col. Temp.:** ambient**Detection:** Mass Spectrometer (MS) @ amu (ambient)**Analyst Note:** SecurityGuardTM Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, POLAR-RP 4 x 2.0mm, 10/Pk Part No.: AJ0-6075

- Holder Part No.: KJ0-4282



15532

**ANALYTES:**

- 1 beta-Zearalenol
- 2 alpha-Zearalenol
- 3 Zearalenone

