

Fumonisin B1 using Gemini**Column:** Gemini® 5 µm C18 110 Å, LC Column 150 x 3 mm, Ea**Dimensions:** 150 x 3 mm ID**Order No:** 00F-4435-Y0**Elution Type:** Isocratic**Eluent A:** 0.1% Formic acid in Water**Eluent B:** Acetonitrile

Gradient	Step No.	Time (min)	Pct A	Pct B
Profile:	1	0	50	50

Flow Rate: 0.5 mL/min**Col. Temp.:** ambient**Detection:** Mass Spectrometer (MS) @ amu (ambient)**Analyst Note:** SecurityGuard™ Guard Cartridge System extends column lifetime.

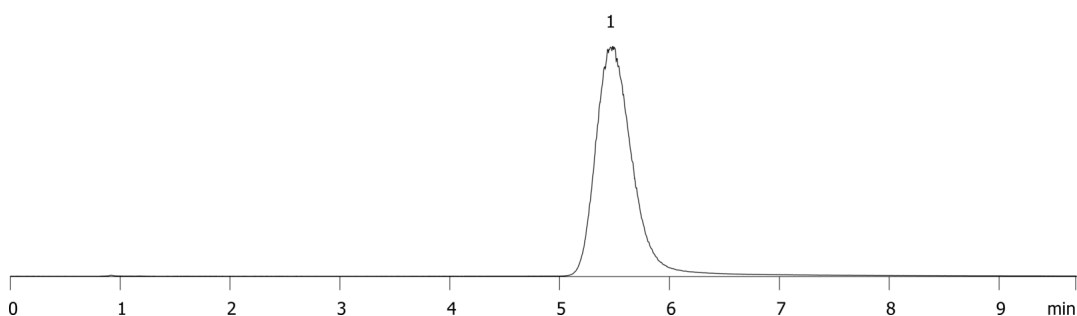
- SecurityGuard Cartridges, Gemini C18 4 x 2.0mm, 10/Pk Part No.: AJ0-7596

- Holder Part No.: KJ0-4282

15537



Products used in this application:

**ANALYTES:****1** Fumonisin B1