

Spearmint Oil on ZB-WAX


Column: Zebtron™ ZB-WAX, GC Cap. Column 20 m x 0.18 mm x 0.18 µm, Ea
Phase: Polyethylene Glycol (PEG)
Dimensions: 20 meters x 0.18 mm x 0.18 µm
Order No: 7FD-G007-08
Oven Profile: 50°C to for 2 min to 100°C at 8°C/min to 225°C at 18°C/min for 4 min

Carrier Gas: Constant Flow Helium, 1 mL/min
Injection: Split 100:1 0.2 µL @ 180°C
Detection: Mass Selective (MSD) (180°C)

15559



Products used in this application:

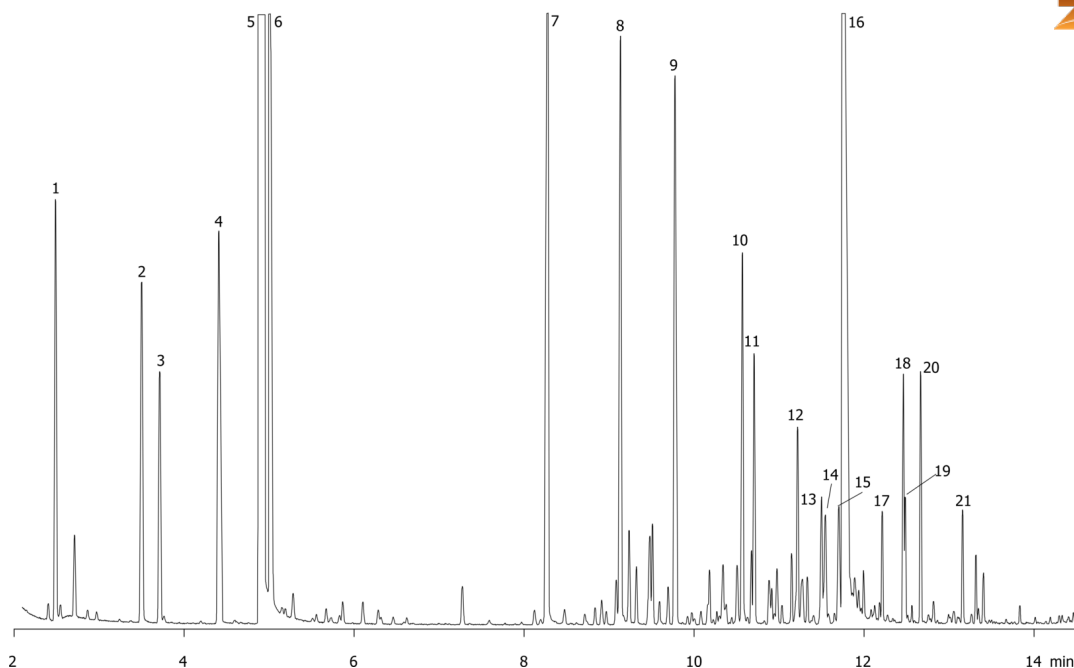

Spearmint Oil on ZB-WAX

ANALYTES:

- 1 alpha-Pinene
- 2 beta-Pinene
- 3 Sabinene
- 4 beta-Myrcene
- 5 D-Limonene
- 6 Eucalyptol
- 7 3-Octanol
- 8 Menthone
- 9 beta-Bourbonene
- 10 Caryophyllene
- 11 Dihydroarvone
- 12 beta-Farnesene
- 13 alpha-Terpineol
- 14 Germacrene D
- 15 Piperitone
- 16 Carvone
- 17 trans-Carvone oxide
- 18 trans-Carveol
- 19 cis-Carvone Oxide
- 20 cis-Carveol
- 21 cis-Jasmone

