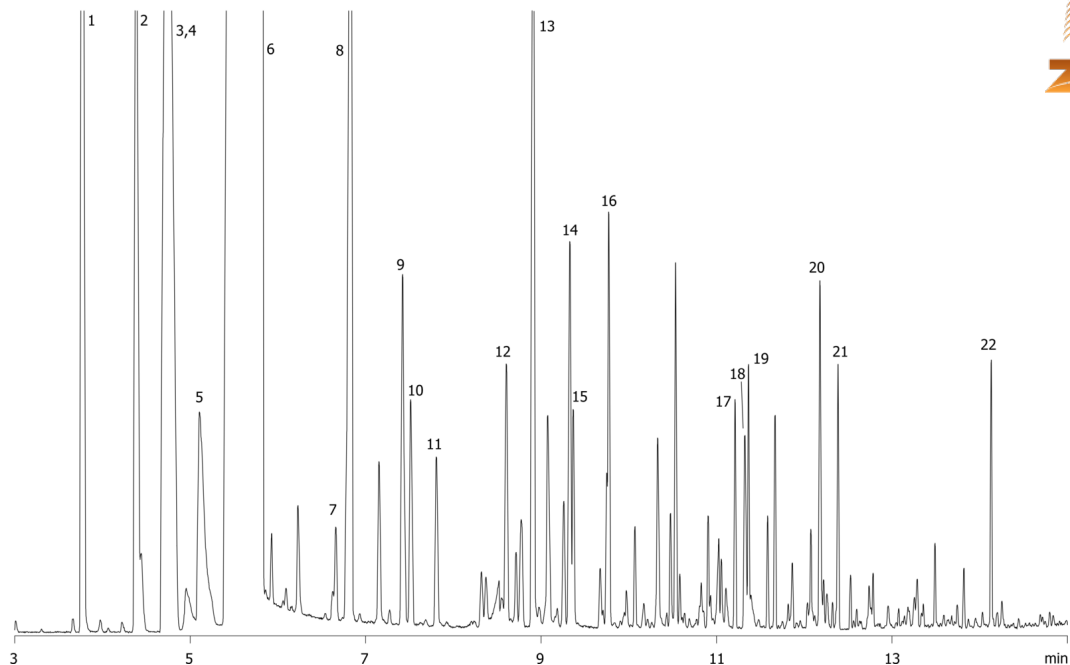


Cold-Pressed Orange Oil on ZB1ms 20x0.18x0.18**Column:** Zebron[™] ZB-1MS, GC Cap. Column 20 m x 0.18 mm x 0.18 µm, Ea**Phase:** 100% Dimethylpolysiloxane**Dimensions:** 20 meters x 0.18 mm x 0.18 µm**Order No:** 7FD-G011-08**Oven Profile:** 50°C to 100°C at 6°C/min to 275°C at 18°C/min for 4 min**Carrier Gas:** Constant Flow Helium, 1 mL/min**Injection:** Split 100:1 1 µL @ 180°C**Detection:** Mass Spectrometer (MS) (240°C)

15571

**Products used in this application:**

Cold-Pressed Orange Oil on ZB1ms 20x0.18x0.18

ANALYTES:

- 1 α -Pinene
- 2 Sabinene
- 3 beta-Myrcene
- 4 Octanal
- 5 3-Carene
- 6 Limonene
- 7 Nonanal
- 8 Linalool
- 9 cis-Limonene oxide
- 10 trans-Limonene oxide
- 11 Citronellal
- 12 α -Terpineol
- 13 Decanal
- 14 Neral
- 15 Carvone
- 16 Geranial
- 17 α -Cubebene
- 18 β -Cubenene
- 19 Dodecanal
- 20 Valencene
- 21 Cadinine
- 22 Nootkatone

