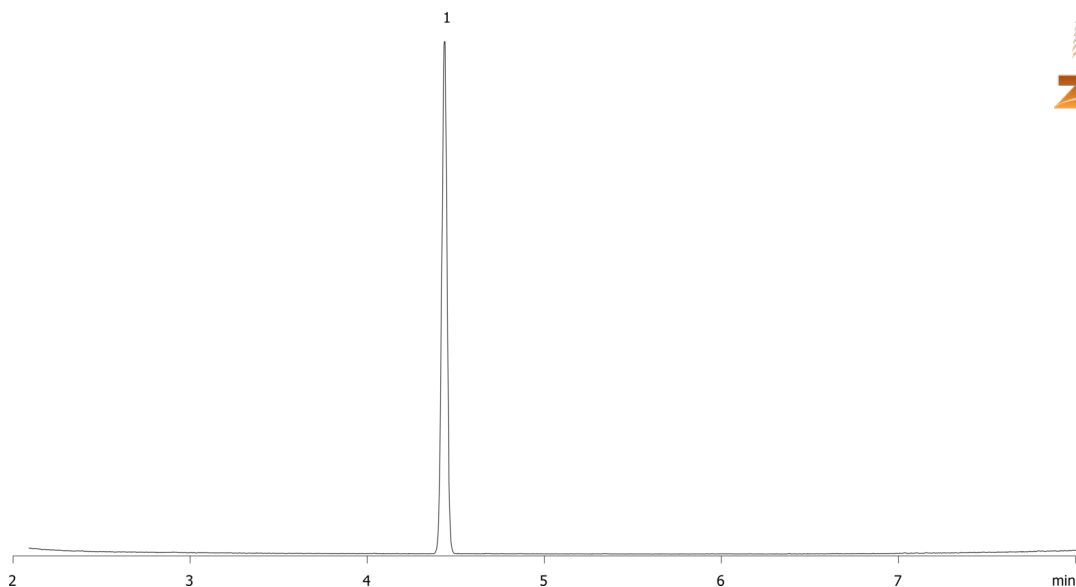


**PCP (Phencyclidine) on ZB1ms 20x0.18x0.18****Column:** Zebtron<sup>™</sup> ZB-1MS, GC Cap. Column 20 m x 0.18 mm x 0.18 µm, Ea**Phase:** 100% Dimethylpolysiloxane**Dimensions:** 20 meters x 0.18 mm x 0.18 µm**Order No:** 7FD-G011-08**Oven Profile:** 175°C to 215°C at 8°C/min to 300°C at 25°C/min**Carrier Gas:** Constant Flow Helium, 0.8 mL/min**Injection:** Split 10:1 1 µL @ 290°C**Detection:** Mass Selective (MSD) (240°C)

15577

**Products used in this application:****ANALYTES:****1** PCP