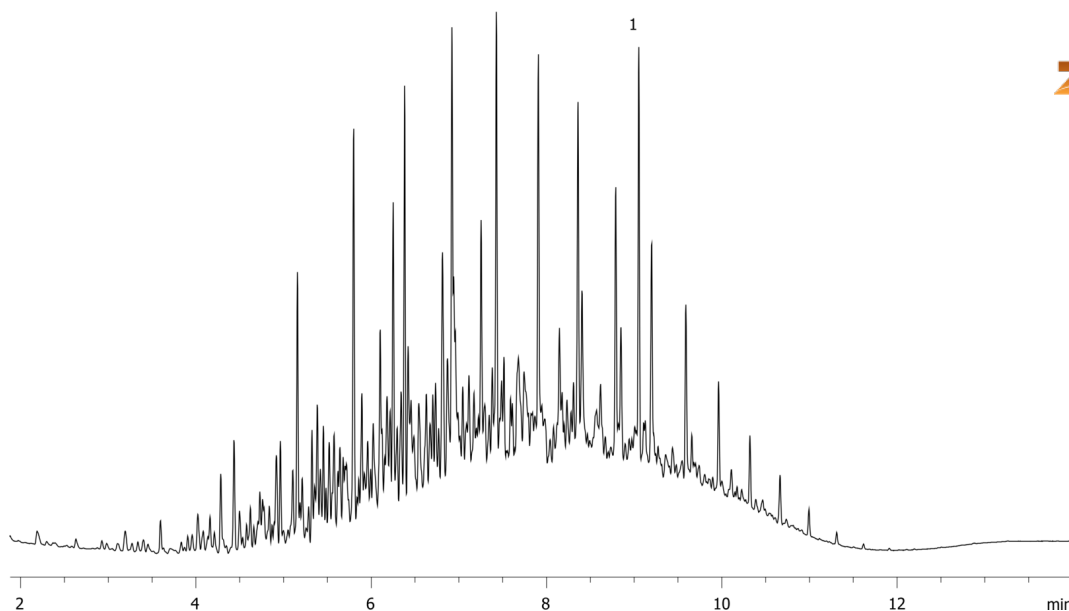


Diesel Fuel Oil #2 on ZB1ms with 5m Guardian ZB1ms 30+5x0.25x0.25**Column:** Zebtron™ ZB-1MS w/GUARDIAN™, 5 m, GC Cap. Column 30 m x 0.25 mm x 0.25 µm, Ea**Phase:** 100% Dimethylpolysiloxane**Dimensions:** 30 meters x 0.25 mm x 0.25 µm**Order No:** 7HG-G011-11-GGA**Oven Profile:** 50°C for 2 min to 320°C at 25°C/min for 5 min**Carrier Gas:** Constant Flow Helium, 4.5 mL/min**Injection:** Split 10:1 1 µL @ 300°C**Detection:** Flame Ionization (FID) (320°C)

15579

**Products used in this application:****ANALYTES:****1** o-Terphenyl