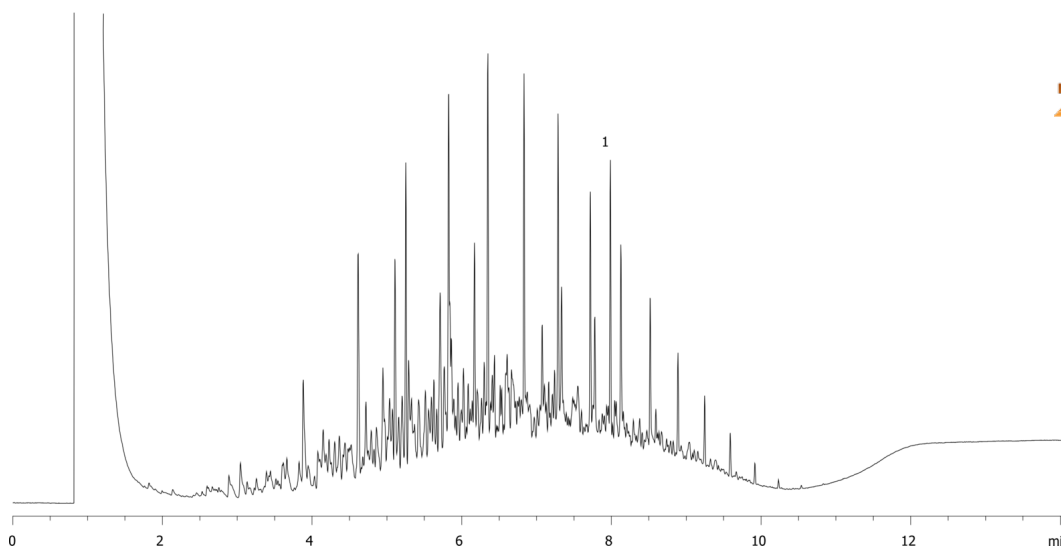


Diesel Fuel Oil #2 on ZB1ms 20x0.18x0.18**Column:** Zebtron[™] ZB-1MS, GC Cap. Column 20 m x 0.18 mm x 0.18 µm, Ea**Phase:** 100% Dimethylpolysiloxane**Dimensions:** 20 meters x 0.18 mm x 0.18 µm**Order No:** 7FD-G011-08**Oven Profile:** 100°C for 2min to 320°C at 25°C/min for 5 min**Carrier Gas:** Constant Flow Helium, 1.2 mL/min**Injection:** Split 10:1 1 µL @ 300°C**Detection:** Flame Ionization (FID) (320°C)

15582

**Products used in this application:****ANALYTES:****1** o-Terphenyl