

Diazepam Analysis in Urine

Column: Onyx™ Monolithic C18, LC Column 100 x 4.6 mm, Ea

Dimensions: 100 x 4.6 mm ID

Order No: CH0-7643

Elution Type: Gradient

Eluent A: 0.1% formic acid in Water

Eluent B: 0.1% formic acid in Acetonitrile

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	95	5
	2	4	10	90

Flow Rate: 4 mL/min


Col. Temp.: 30 °C

Detection: Mass Spectrometer (MS) @ (ambient)

15588



Products used in this application:

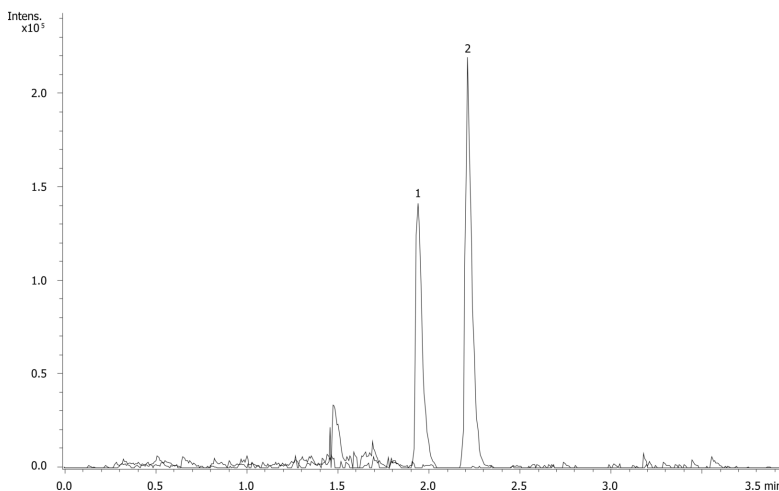

ANALYTES:

- 1 Nordiazepam
- 2 Diazepam

