

'Soy Complex' Extract (Finished Product)

Column: Onyx™ Monolithic C18, LC Column 100 x 4.6 mm, Ea

Dimensions: 100 x 4.6 mm ID

Order No: CH0-7643

Elution Type: Gradient

Eluent A: 0.1% TFA in Water

Eluent B: 0.08% TFA in Acetonitrile

Gradient	Step No.	Time (min)	Pct A	Pct B
Profile:	1	0	95	5
	2	15	30	70

Flow Rate: 1 mL/min


Col. Temp.: 30 °C

Detection: UV-Vis Abs.-Variable Wave.(UV) @ 280 nm (ambient)

Analyst Note: 2 columns coupled end to end
15593



Products used in this application:

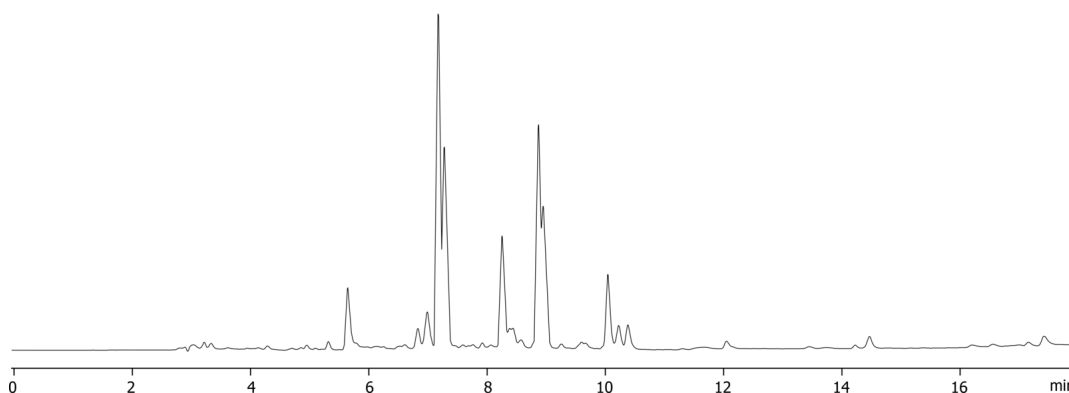

ANALYTES:

- 1 Soy complex

