

Ibuprofen Analysis in Urine

Column: Onyx™ Monolithic C18, LC Column 100 x 4.6 mm, Ea

Dimensions: 100 x 4.6 mm ID

Order No: CH0-7643

Elution Type: Gradient

Eluent A: 0.1% Formic acid in Water

Eluent B: 0.1% Formic acid in Acetonitrile

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	95	5
	2	4	10	90

Flow Rate: 4 mL/min


Col. Temp.: 30 °C

Detection: Mass Spectrometer (MS) @ _ (ambient)

15594



Products used in this application:

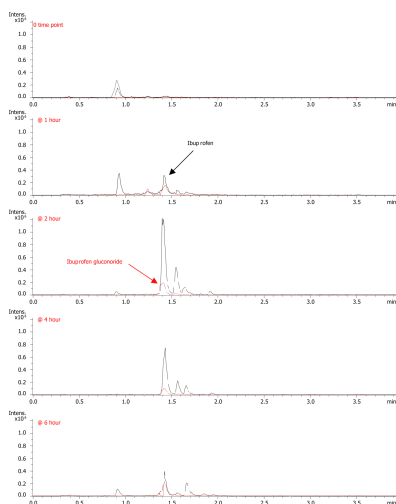

ANALYTES:

1 Ibuprofen

