

Mycotoxin - Desoxynivalenol and Zearalenon**Column:** SynergiTM 4 μ m Polar-RP 80 Å, LC Column 150 x 2 mm, Ea**Dimensions:** 150 x 2 mm ID**Order No:** 00F-4336-B0**Elution Type:** Gradient**Eluent A:** Methanol**Eluent B:** Water with 5% Methanol

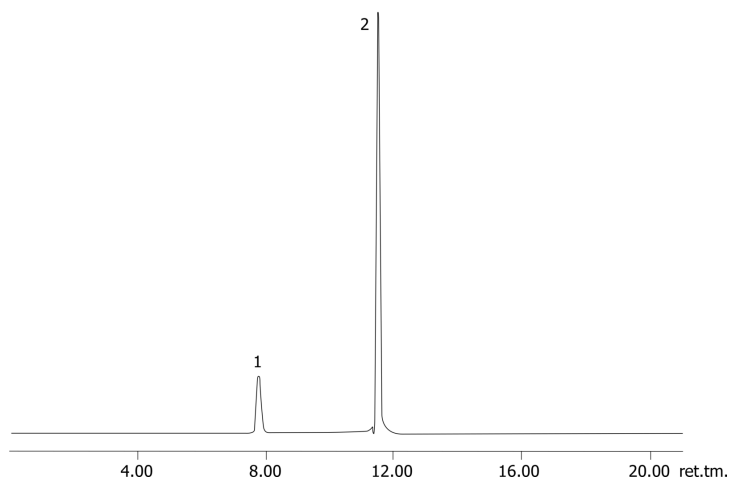
Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	10	90
	2	5	90	10
	3	10	90	10
	4	15	10	90
	5	20	10	90

Flow Rate: 0.2 mL/min**Col. Temp.:** 45 °C**Detection:** Mass Spectrometer (MS) @ amu (ambient)**Detector Info:** MS/MS Kodiak 1500**Analyst Note:** SecurityGuardTM Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, POLAR-RP 4 x 2.0mm, 10/Pk Part No.: AJ0-6075

- Holder Part No.: KJ0-4282

15631

**ANALYTES:**

- 1 Desoxynivalenol
- 2 Zearalenone

