

Performance test on Onyx C18 (column #3)

Column: Onyx™ Monolithic C18, LC Column 100 x 4.6 mm, Ea

Dimensions: 100 x 4.6 mm ID

Order No: CH0-7643

Elution Type: Isocratic

Eluent A: Water

Eluent B: Acetonitrile

Gradient

Step No.	Time (min)	Pct A	Pct B
1	0	40	60

Profile:


Flow Rate: 2 mL/min

Col. Temp.: ambient

Detection: UV-Vis Abs.-Variable Wave.(UV) @ 254 nm (ambient)
15644



Products used in this application:

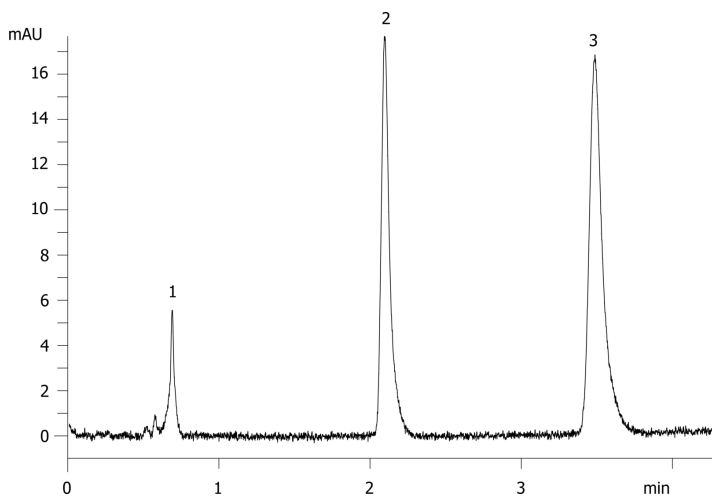

ANALYTES:

1 Thiourea

Retention Time: 0.69 min

2 Progesterone

Retention Time: 2.097 min

3 Anthracene

Retention Time: 3.485 min

