

Performance test on Onyx C18 (short column)


Column: Onyx™ Monolithic C18, LC Column 50 x 4.6 mm, Ea
Dimensions: 50 x 4.6 mm ID
Order No: CH0-7644
Elution Type: Isocratic
Eluent A: water
Eluent B: acetonitrile
Gradient Profile:

Step No.	Time (min)	Pct A	Pct B
1	0	40	60

Flow Rate: 2 mL/min
Col. Temp.: ambient
Detection: UV-Vis Abs.-Variable Wave.(UV) @ 254 nm (ambient)



Products used in this application:

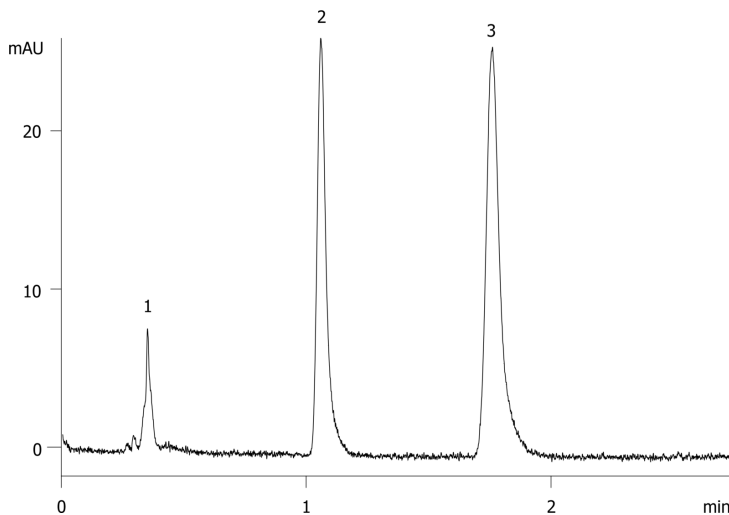

ANALYTES:

- 1** Thiourea
Retention Time: 0.353 min
- 2** Progesterone
Retention Time: 1.06 min
- 3** Anthracene
Retention Time: 1.759 min

